



TECNILABO

HS Code: 84184080
KUB4075 PLL ADV
LT Freezer -45°C

Technical Features

Net Storage Volume	86 lt (3,04 cu ft)
2" Cryoboxes Max Capacity	48
1,5/2ml Vials Max Capacity	4.800
Racks Max Capacity	6
Temperature Control Range	-20°C/-45°C
Factory T Set Point	-40°C
Ambient Operating Temperature	+10°C/+35°C, 50% RH (N-ST Climate Class)
Noise Level	< 42 dBA

Power Supply	230V ±10% / 50Hz
Standard Electrical Plug	Industrial plug CEE 16A, 2P+E
Max Power Absorption	4A
Daily Energy Consumption*	3,80 kWh/24h



External Dimensions	W 990 x D 635 x H 780 mm W 38,97 x D 25,00 x H 30,71 inch.
Internal Dimensions	W 428 x D 343 x H 588 mm W 16,85 x D 13,50 x H 23,15 inch.
Weight	75 kg

Packing Dimensions	Carton Box + Pallet: W 1100 x D 800 x H 1000 mm, 85 kg Wooden Crate: W 1200 x D 900 x H 1100 mm, 125 kg
--------------------	--

External Surface	Anti-corrosion hot-dip galvanized steel sheet, covered with PVC film
Internal Surface	AISI 304 18/10 stainless steel
External/Internal Edges	Rounded to facilitate cleaning
Insulation Type	High-density CFC-free PUR foam
Insulation Thickness	90 mm on sides / 90 mm rear / 80 mm outer door
Door Type	Solid, with n°2 heavy-duty side hinges for the max stability
Door Sealing	Triple-layer silicone gasket, heated by the refrigerant to avoid the formation of ice
Handle	Ergonomic with integrated key-lock system (n°2 keys included)


Standard Configuration	n°1 shelf in stainless steel (AISI 304), adjustable in height (adjustment steps of 1.5 cm) with max load of 55kg each n°1 insulated inner door (inner chamber divided into 2 compartments); inner door can be removed for cleaning n°4 leveling feet Porthole (Ø 25 mm) for introducing in the chamber independent probes connected to external monitoring system/BMS External vacuum relief valve (automatic & manual pressure release); the valve can be removed for cleaning operations Back-up battery with self-charging circuit to ensure recording of T and alarm data in case of power failure (48 h autonomy) Volt-free contact for connection to an external alarm monitoring system and BMS
------------------------	--

Cooling System

Cooling	Static cooling with n°1 hermetic compressor (with thermal protector); air-cooled condenser with filterless design
HC Refrigerant (CFC/HCFC free)	R452a
Defrost System	Manual

PL ADV Control System

Controller	Microprocessor Controller
Display	LED display with 5-keys capacitive keyboard, with 0.1°C resolution
Temperature Probes	n°2 NTC probes (for evaporator / temperature control) + n°1 Pt1000 (for monitoring)
Access Control	3-levels password security system (User/Service/Admin)
Connectivity	In-built Wi-Fi module for remote connection to mobile APP ConneX (for IOS & Android)/Cloud The main screen displays the temperature set point, functional icons, current time, alarm events Data logger function with automatic recording of temperature data (at 30-seconds intervals) and events
Features	Flash memory capable of recording data for a minimum period of 1 year In-built power stabilizer with linear voltage regulation (± 20V) Possibility of remotely setting and modifying the device parameters and to view the temperature graph in real time Possibility to download data in CSV/PDF format via APP on smartphone/tablet/PC (no need for USB) Automatic sending of push notifications (via APP) and/or email (via Cloud) in the event of an alarm Possibility of setting automatic sending of weekly temperature and event reports via email
Remote Management and Monitoring with IOS/Android ConneX APP and/or KW Cloud	- Power failure - Door open (interval can be freely set by the user) - High/Low T (limits can be set freely by the user) - Temperature probes failure - Back-up battery failure / low charge - Dirty condenser / High ambient temperature
Main Alarms List (Audio/visual notification)	

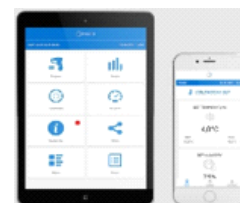

Accessories Available On Request


Frontal extraction racks made in AISI 304 stainless steel for cryoboxes storage
 n°4 swivel castors (n°2 front wheels with brake)
 Additional stainless steel (AISI 304) shelf
 Voltage stabilizer
 Electrical plug different from standard
 Packing in wooden crate



PL ADV controller



Remote alarm contacts

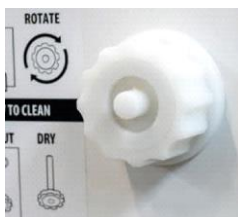


Vacuum relief valve



Ergonomic handle with key-lock system



Shelves adjustment system

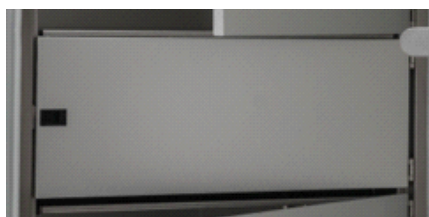

Insulated inner door



Heated triple-layer silicone gasket

* with +20°C ambient temperature



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.