




Overview: model class **Q-CLASSIC-30**

Safe and approved storage of flammable hazardous substances in working areas



Price-performance winner

- the sophisticated design and the use of melamine-resin-coated, fireproof special plates ensures an optimum price/performance ratio for the Q-CLASSIC-30

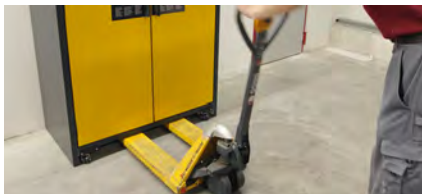

High flexibility

- integrated transport base for simple inhouse transport
- cabinet can be transported upright through doors (all Q-LINE models)
- adjustable feet installed in the base allow the safe and fast compensation of minor surface unevenness



Warranty up to 5 years


- quality from which you benefit
- you get an extended manufacturer's warranty of up to 5 years on the Q-CLASSIC-30 cabinets

Further features of this model class



Modern bicolour design

- the Q-CLASSIC brings colour to your workplace: choice of seven door colours at no extra cost
- language-free labelling for optimum international usability

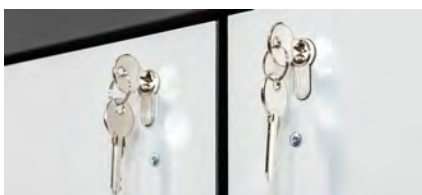

Cylinder lock

- door lockable with cylinder lock, the cabinets can be integrated by the customer in an existing locking system

Possible interior equipment



- metal-free interior equipment** for the storage of aggressive hazardous materials
- as standard with **removable plastic sump**



- drawers** with high storage volumes, ideal for storing small containers



- height-adjustable shelves** in different versions



- bottom collecting sump** – safe collection of leakages
- usable as storage level in conjunction with perforated insert



Safety storage underbench cabinets for the storage of flammable hazardous materials in working areas according to EN 14470-1 and TRGS 510



CLASSIC



- **fire resistance of 30 minutes** (type 30), type-tested in accordance with EN 14470-1



- **GS-tested** (EN 14470-1, EN 14727, document EK5/AK4 09-10, ProdSG)
- **50,000** wear-free opening and closing actions
- test basis: all cabinets of this model group correspond to the requirements of the Equipment and Product Safety Act (§7 Para. 1) with respect to the guarantee of safety and health



- all cabinets of this model group have been successfully tested in accordance with **FM test standard 6050** and are compliant to NFPA Code 30, NFPA 1 Fire Code, OSHA



- **excellent products with special manufacturer warranty**
- extended **manufacturer warranty of up to 5 years** in combination with the asecos BASICPlus tariff





Q-CLASSIC-30 ◀ 116 cm ▶




Safety storage cabinet

Q-CLASSIC-30 model Q30.195.116

Body colour Anthracite grey, with wing doors in Warning yellow (RAL 1004), interior equipment with 3 x shelf, 1 x perforated insert, 1 x bottom collecting sump (sheet steel powder-coated)

Order No. 31003-060-33635

(base cover and recirculating air filter system optional)

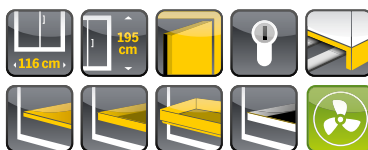

Safe and approved storage of flammable hazardous substances in working areas

Function / construction:

- **robust construction and longevity:** triple hinged door, cabinet body made of melamine resin and flame retardant special plates, safety elements assembled outside the storage compartment for increased protection against corrosion, doors with metal front as well as scratch- and impact-resistant surface, easy to clean
- **easy handling with comfort:** doors stop in any position, self-closing in the event of fire
- **no unauthorised use:** doors lockable with profile cylinder (integration in an existing locking system possible)
- **easy transport:** integrated transport base with removable, optional base cover or adjusting aids with base cover

- **ventilation:** integrated air ducts ready for connection (DN 75) to a technical exhaust system

Available equipment

- shelves with load capacity of 75 kg
- tray shelves with load capacity of 150 kg
- perforated insert useable as storage area
- bottom collecting sump
- drawers – high storage capacity, optimal use of the available storage capacities
- drawers heavy-duty with load capacity of 100 kg
- metal-free tray shelves and collecting sumps
- storage box
- sump inliner (PE or PP)



Maximum flexibility

No problem with the Q-CLASSIC-30 including transport base and technical ventilation with recirculating air filter system (optional), Q-Mover (optional) for in-house transport



Drawers

High storage capacity
Ideal for the storage of small containers



SUCCESSFULLY
TESTED IN
ACCORDANCE
WITH FM TEST
STANDARD 6050



Extraction unit

Safe ventilation of safety storage cabinets

plug-in – fast and easy to install

Further information and products from page 397

Order No. 14218



Recirculating air filter system

Avoid complex wall break-throughs and expensive extraction piping

Further information from page 393

Order No. 25099

**Safety storage cabinet****Q-CLASSIC-30 model Q30.195.116**

Body colour Anthracite grey, with wing doors in Light grey (RAL 7035), interior equipment with 4 x drawer (sheet steel powder-coated)

Order No. 31003-061-33647

(base cover optional)

**Safety storage cabinet****Q-CLASSIC-30 model Q30.195.116**

Body colour Anthracite grey, with wing doors in Gentian blue (RAL 5010), interior equipment with 6 x drawer (sheet steel powder-coated)

Order No. 31003-067-33649

(base cover and recirculating air filter system optional)

7 door colours are available **at no extra cost:**



Order example - how to order correctly



	Model	Information / equipment	Order No.
	Q30.195.116	without interior equipment	31003

	Body colour	Door colour	Order No.
	Anthracite grey	Warning yellow RAL 1004	060
	Anthracite grey	Light grey RAL 7035	061
	Anthracite grey	Traffic red RAL 3020	066
	Anthracite grey	Gentian blue RAL 5010	067
	Anthracite grey	Reseda green RAL 6011	068
	Anthracite grey	asecos silver	069
	Anthracite grey	Pure white RAL 9010	070



	Interior equipment packages (can only be ordered in combination with the cabinet)	Material	Order No.
			old code new code
	3x shelf, 1x perforated insert, 1x bottom collecting sump (V=33.0L)	sheet steel powder-coated RAL 7035	30003 33635
	4x shelf, 1x perforated insert, 1x bottom collecting sump (V=33.0L)	sheet steel powder-coated RAL 7035	30004 33636
	3x shelf, 1x bottom collecting sump (V=33.0L)	sheet steel powder-coated RAL 7035	30005 33637
	4x shelf, 1x bottom collecting sump (V=33.0L)	sheet steel powder-coated RAL 7035	30006 33638
	3x tray shelf (V=33.0L), 1x perforated insert, 1x bottom collecting sump (V=55.5L)	sheet steel powder-coated RAL 7035	30007 33640
	3x shelf, 1x perforated insert, 1x bottom collecting sump (V=33.0L)	stainless steel	30008 33642
	4x shelf, 1x perforated insert, 1x bottom collecting sump (V=33.0L)	stainless steel	30009 33643
	3x shelf, 1x bottom collecting sump (V=33.0L)	stainless steel	30010 33644
	4x shelf, 1x bottom collecting sump (V=33.0L)	stainless steel	30011 33645
	3x tray shelf (V=33.0L), 1x bottom collecting sump (V=33.0L)	PP grey	HIGHLIGHT 30013 30013
	4x tray shelf (V=33.0L), 1x bottom collecting sump (V=33.0L)	PP grey	HIGHLIGHT 30014 30014
	4x drawer (V=30.0L)	sheet steel powder-coated RAL 7035	30015 33647
	5x drawer (V=30.0L)	sheet steel powder-coated RAL 7035	30016 33648
	6x drawer (V=30.0L)	sheet steel powder-coated RAL 7035	30017 33649
	4x drawer (V=30.0L)	stainless steel	30018 33651
	5x drawer (V=30.0L)	stainless steel	30019 33652
	6x drawer (V=30.0L)	stainless steel	30020 33653



Individual configuration of model Q30.195.116

Model	Information / equipment	W x D x H (mm)	Order No.
Q30.195.116	without interior equipment	1164 x 620 x 1947	31003


Body colour	Door colour	Order No.
Anthracite grey	Warning yellow RAL 1004	060 4
Anthracite grey	Light grey RAL 7035	061 4
Anthracite grey	Traffic red RAL 3020	066 4
Anthracite grey	Gentian blue RAL 5010	067 4
Anthracite grey	Reseda green RAL 6011	068 4
Anthracite grey	asecos silver	069 4
Anthracite grey	Pure white RAL 9010	070 4

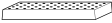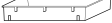Interior equipment	Material	W x D x H (mm)	 1)	 2)	 3)	Order No.
	shelf	sheet steel powder-coated RAL 7035	1040 x 432 x 30*	75.00		old code 30021 new code 30021 2
		stainless steel	1038 x 425 x 29*	75.00		30024 30024 2
	tray shelf	sheet steel powder-coated RAL 7035	1025 x 445 x 71*	150.00	33.00 30.00	30028 30028 2
		stainless steel	1025 x 445 x 71*	150.00	33.00 30.00	30064 30064 2
	tray shelf, metal-free	PP grey	1002 x 472 x 74*	50.00	33.00 30.00	30027 30027 2
	perforated insert	sheet steel powder-coated RAL 7035	1024 x 417 x 60	75.00		5536 5536 2
		stainless steel	1024 x 417 x 60	75.00		5773 5773 2
	perforated insert, height = 135 mm	sheet steel powder-coated RAL 7035	1024 x 417 x 135	135.00		14084 14084 2
	bottom collecting sump	sheet steel powder-coated RAL 7035	1024 x 489 x 86*		33.00 30.00	30022 33634 2
		stainless steel	1024 x 489 x 86*	150.00	33.00 30.00	30025 33641 2
	bottom collecting sump, V = 55,5 ltr.	sheet steel powder-coated RAL 7035	1017 x 482 x 127*		55.50 50.45	30023 33639 2
		stainless steel	1018 x 438 x 120*		55.50 50.00	30026 30026 10
	bottom collecting sump, metal-free	PP grey	1002 x 472 x 74*		33.00 30.00	26989 26989 2
	drawer	sheet steel powder-coated RAL 7035	844 x 488 x 90*	60.00	30.00 27.27	30013 33646 2
		stainless steel	844 x 488 x 90*	60.00	30.00 27.27	30033 33650 2
	drawer (heavy-duty)	sheet steel powder-coated RAL 7035	843 x 487 x 162*	100.00	66.00	30032 30032 4
		stainless steel	843 x 487 x 162*	100.00	66.00	30034 30034 10
	storage box, left hinged	sheet steel powder-coated RAL 7035	380 x 420 x 360			30029 30029 2
	storage box, right hinged	sheet steel powder-coated RAL 7035	380 x 420 x 360			30030 30030 2
	sump inliner, for shelf	polyethylene black	1000 x 392 x 50*			7223 7223 2
	sump inliner, for drawer	PP grey	834 x 450 x 96*			10054 10054 2
	anti-slip rubber mat, for drawer	black	840 x 490 x 1			8699 8699 2
	pipe lead-through	sheet steel powder-coated RAL 7035	182 x 104 x 48			30035 30035 2
	base cover, with adjusting aids	panel melamine resin-coated RAL 7016				31005 31005 2
	base cover	panel melamine resin-coated RAL 7016				31925 31925 2
	Q-mover, 2 pieces	steel powder-coated textured RAL 7016	600 x 140 x 130	450.00		30037 30037 2
	box for chemicals	PP blue	410 x 300 x 260		18.00 16.36	30038 30038 2

Accessories	Order No.
recirculating air filter system with exhaust air monitoring	25099 2
extraction unit without exhaust air monitoring	14218 2

* usable storage space W x D [mm]


1) load capacity [kg], 2) retention capacity [ltr], 3) maximum container volume in the cabinet [ltr]

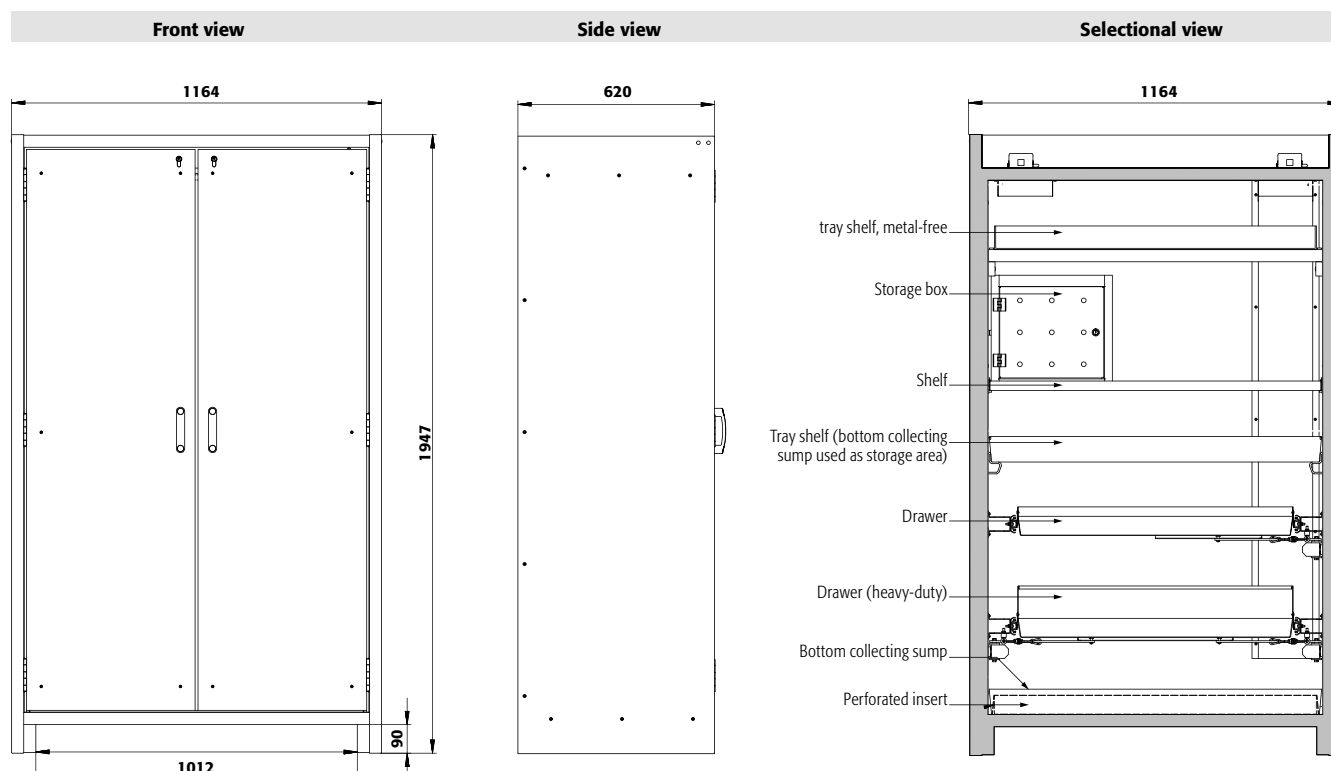**Technical data**

External dimensions W x D x H	mm	1164 x 620 x 1947
Internal dimensions W x D x H	mm	1050 x 519 x 1626
Weight without interior equip.	kg	240
Air change 10 times	m ³ /h	9
Diff. pressure (air change 10 times)	Pa	10
Diff. pressure (air change 10 times EK5/AK4)	Pa	10
Maximum load	kg	600
Distributed load	kg/m ²	445

Transport base

Entry height transport base	mm	90
Entry width transport base	mm	1012

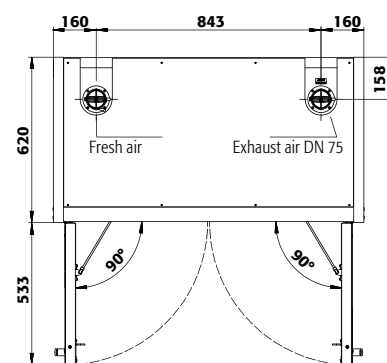Q30.195.116**Top view****Q30.195.116****Drawers/number**


	9	8	7	6	5	4+1	4	3	2
Distance ⑨ (mm)	130								
Distance ⑧ (mm)	135	130							
Distance ⑦ (mm)	135	170	160						
Distance ⑥ (mm)	135	170	200	225					
Distance ⑤ (mm)	135	170	200	230	225	225			
Distance ④ (mm)	135	170	200	230	295	230	320		
Distance ③ (mm)	135	170	200	230	295	230	360	420	
Distance ② (mm)	170	170	200	230	295	230	360	520	865
Distance ① (mm)	200	200	200	230	295	490	390	520	615

Q30.195.116**Shelves/number**

	4	3
Distance ④ (mm)	370	
Distance ③ (mm)	285	400
Distance ② (mm)	315	380
Distance ① (mm)	315	410

Standard storage heights of the shelves**Q30.195.116****Quantity**

Quantity	Type of container	Container size
90	Aerosol cans	400 ml
40	Cans	1 litre
40	Laboratory glass bottles	1 litre
40	Laboratory PE-bottles	1 litre
15	Laboratory glass bottles	2.5 litres

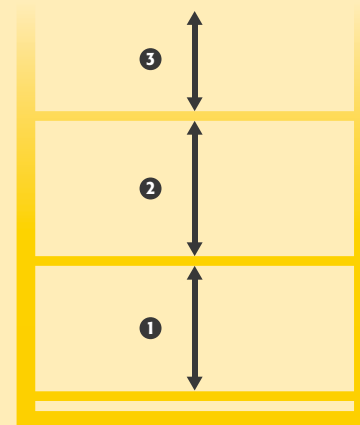Storage capacities per drawer/type of container



Q-CLASSIC-30 ◀ 086 cm ▶




Safety storage cabinet

Q-CLASSIC-30 model Q30.195.086.WD

Body colour Anthracite grey, with wing door in Warning yellow (RAL 1004), interior equipment with 3 x shelf, 1 x perforated insert, 1 x bottom collecting sump (sheet steel powder-coated)

Order No. 31008-060-33622 4

(base cover optional)

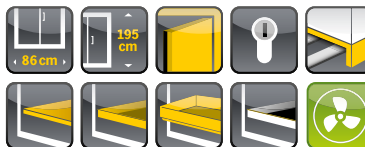

Safe and approved storage of flammable hazardous substances in working areas

Function / construction:

- **robust construction and longevity:** triple hinged door, cabinet body made of melamine resin and flame retardant special plates, safety elements assembled outside the storage compartment for increased protection against corrosion, doors with metal front as well as scratch-and impact-resistant surface, easy to clean
- **easy handling with comfort:** door stops in any position, self-closing in the event of fire
- **no unauthorised use:** door lockable with profile cylinder (integration in an existing locking system possible)
- **easy transport:** integrated transport base with removable, optional base cover or adjusting aids with base cover

- **ventilation:** integrated air ducts ready for connection (DN 75) to a technical exhaust system

Available equipment

- shelves with load capacity of 75 kg
- tray shelves with load capacity of 150 kg
- perforated insert useable as storage area
- bottom collecting sump
- drawers – high storage capacity, optimal use of the available storage capacities
- drawers heavy-duty with load capacity of 100 kg
- metal-free tray shelves and collecting sumps
- storage box
- sump inliner (PE or PP)



Drawers

High storage capacity
Ideal for the storage of small containers



Metal-free interior equipment

Ideal for the storage of flammable, aggressive hazardous substances. Tray shelves with chemical resistant surface lamination and removeable plastic sump, bottom collecting plastic sump.



**LARGE STORAGE
CAPACITY**
—
**OUTSTANDING PRICE
PERFORMANCE**



Shelf

Height adjustable (32 mm increments), with 30 mm upstand



**SUCCESSFULLY
TESTED IN
ACCORDANCE
WITH FM TEST
STANDARD 6050**

**Safety storage cabinet****Q-CLASSIC-30 model Q30.195.086.WD**

Body colour Anthracite grey, with wing door in Traffic red (RAL 3020), interior equipment with 4 x drawer (sheet steel powder-coated)

Order No. 31008-066-33627

(base cover and extraction unit optional)

Safety storage cabinet**Q-CLASSIC-30 model Q30.195.086.WD**

Body colour Anthracite grey, with wing door in Light grey (RAL 7035), interior equipment with 6 x drawer (sheet steel powder-coated)

Order No. 31008-061-33629

(base cover optional)

7 door colours are available **at no extra cost:**



Order example - how to order correctly



Model	Information / equipment	Order No.
Q30.195.086.WD	without interior equipment, left hinged	31008
Q30.195.086.WDR	without interior equipment, right hinged	31009



Body colour	Door colour	Order No.
Anthracite grey	Warning yellow RAL 1004	060
Anthracite grey	Light grey RAL 7035	061
Anthracite grey	Traffic red RAL 3020	066
Anthracite grey	Gentian blue RAL 5010	067
Anthracite grey	Reseda green RAL 6011	068
Anthracite grey	asecos silver	069
Anthracite grey	Pure white RAL 9010	070



Interior equipment packages (can only be ordered in combination with the cabinet)	Material	Order No.
		old code new code
3x shelf, 1x perforated insert, 1x bottom collecting sump (V=22.0L)	sheet steel powder-coated RAL 7035	30043 33622
4x shelf, 1x perforated insert, 1x bottom collecting sump (V=22.0L)	sheet steel powder-coated RAL 7035	30044 33623
3x shelf, 1x bottom collecting sump (V=22.0L)	sheet steel powder-coated RAL 7035	30045 33624
4x shelf, 1x bottom collecting sump (V=22.0L)	sheet steel powder-coated RAL 7035	30046 33625
3x shelf, 1x perforated insert, 1x bottom collecting sump (V=22.0L)	stainless steel	30047 30047
4x shelf, 1x perforated insert, 1x bottom collecting sump (V=22.0L)	stainless steel	30048 30048
3x shelf, 1x bottom collecting sump (V=22.0L)	stainless steel	30049 30049
4x shelf, 1x bottom collecting sump (V=22.0L)	stainless steel	30050 30050
3x tray shelf (V=22.0L), 1x bottom collecting sump (V=22.0L)	PP grey	HIGHLIGHT 30051 30051
4x tray shelf (V=22.0L), 1x bottom collecting sump (V=22.0L)	PP grey	HIGHLIGHT 30052 30052
4x drawer (V=19.0L)	sheet steel powder-coated RAL 7035	30053 33627
5x drawer (V=19.0L)	sheet steel powder-coated RAL 7035	30054 33628
6x drawer (V=19.0L)	sheet steel powder-coated RAL 7035	30055 33629
4x drawer (V=19.0L)	stainless steel	30056 33631
5x drawer (V=19.0L)	stainless steel	30057 33632
6x drawer (V=19.0L)	stainless steel	30058 33633



Individual configuration of model Q30.195.086.WD

 Model	Information / equipment	W x D x H (mm)	Order No.
Q30.195.086.WD	without interior equipment, left hinged	864 x 620 x 1947	31008
Q30.195.086.WDR	without interior equipment, right hinged	864 x 620 x 1947	31009

 Body colour	Door colour	Order No.
Anthracite grey	Warning yellow RAL 1004	060 
Anthracite grey	Light grey RAL 7035	061 
Anthracite grey	Traffic red RAL 3020	066 
Anthracite grey	Gentian blue RAL 5010	067 
Anthracite grey	Reseda green RAL 6011	068 
Anthracite grey	asecos silver	069 
Anthracite grey	Pure white RAL 9010	070 


 Interior equipment	Material	W x D x H (mm)	 1)	 2)	 3)	Order No.	
						old code	new code
 shelf	sheet steel powder-coated RAL 7035	740 x 432 x 30*	75.00			30059	30059 
	stainless steel	740 x 432 x 30*	75.00			30061	30061 
 tray shelf, metal-free	PP grey	702 x 472 x 74*	50.00	23.00	20.91	30063	30063 
 perforated insert	sheet steel powder-coated RAL 7035	742 x 417 x 70	75.00			11165	11165 
	stainless steel	742 x 417 x 70	75.00			11168	11168 
 bottom collecting sump	sheet steel powder-coated RAL 7035	724 x 489 x 83*		22.00	20.00	30060	33621 
 drawer	sheet steel powder-coated RAL 7035	544 x 488 x 90*	60.00	19.00	1727	30065	33626 
 bottom collecting sump	stainless steel	724 x 489 x 83*		23.00	20.91	30062	30062 
 bottom collecting sump, metal-free	PP grey	702 x 472 x 74*		23.00	20.91	26990	26990 
 drawer	stainless steel	544 x 488 x 90*	60.00	19.00	1727	30067	33630 
 drawer (heavy-duty)	sheet steel powder-coated RAL 7035	543 x 487 x 162*	100.00	45.00		30066	30066 
	stainless steel	543 x 487 x 162*	100.00	45.00		30068	30068 
 storage box, left hinged	sheet steel powder-coated RAL 7035	380 x 420 x 360				30029	30029 
 storage box, right hinged	sheet steel powder-coated RAL 7035	380 x 420 x 360				30030	30030 
 sump inliner, for drawer	polyethylene white	536 x 480 x 90*		22.00	20.00	27948	27948 
 sump inliner, for shelf	PP grey	653 x 406 x 88*				9707	9707 
 anti-slip rubber mat, for drawer	black	540 x 490 x 1				15653	15653 
 pipe lead-through	sheet steel powder-coated RAL 7035	182 x 104 x 48				30035	30035 
 base cover, with adjusting device (available for order only in combination with cabinet)	panel melamine resin-coated RAL 7016					31010	31010 
 base cover	panel melamine resin-coated RAL 7016					31924	31924 
 Q-mover, 2 pieces	steel powder-coated textured RAL 7016	600 x 140 x 130	450.00			30037	30037 
 box for chemicals	PP blue	410 x 300 x 260		18.00	16.36	30038	30038 

 Accessories	Order No.
recirculating air filter system with exhaust air monitoring	25099 
extraction unit without exhaust air monitoring	14218 

* usable storage space W x D (mm)

1) load capacity (kg), 2) retention capacity (ltr), 3) maximum container volume in the cabinet (ltr)

CERTIFIED
SAFETY 
YOU CAN RELY ON.

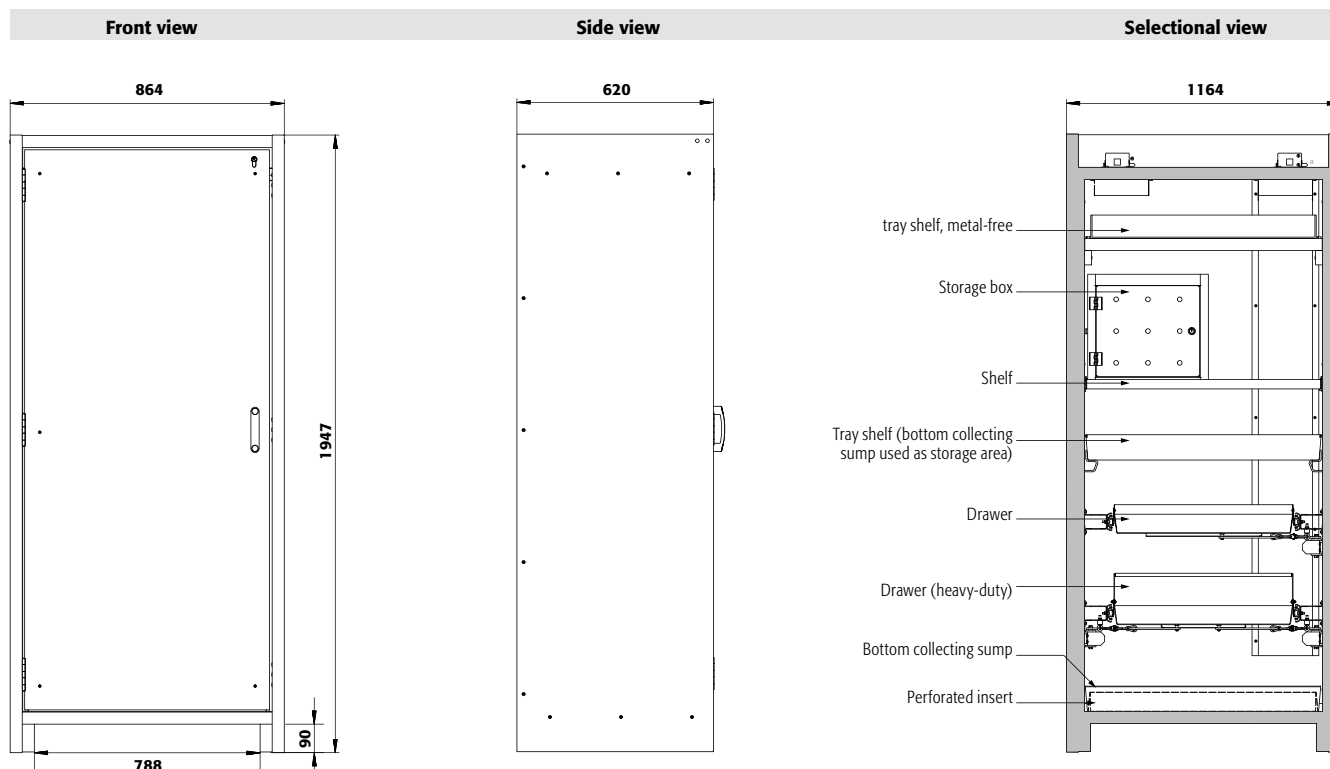

Technical data


Q30.195.086.WD

External dimensions W x D x H	mm	864 x 620 x 1947
Internal dimensions W x D x H	mm	750 x 519 x 1626
Weight without interior equip.	kg	195
Air change 10 times	m ³ /h	6
Diff. pressure (air change 10 times)	Pa	10
Diff. pressure (air change 10 times EK5/AK4)	Pa	10
Maximum load	kg	600
Distributed load	kg/m ²	568

Transport base

Entry height transport base	mm	90
Entry width transport base	mm	788

Top view

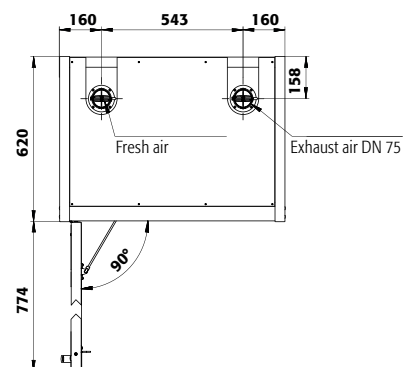

Q30.195.086.WD

Drawers/number

	9	8	7	6	5	4+1	4	3	2
Distance ⑨ (mm)	130								
Distance ⑧ (mm)	135	130							
Distance ⑦ (mm)	135	170	160						
Distance ⑥ (mm)	135	170	200	225					
Distance ⑤ (mm)	135	170	200	230	225	225			
Distance ④ (mm)	135	170	200	230	295	230	320		
Distance ③ (mm)	135	170	200	230	295	230	360	420	
Distance ② (mm)	170	170	200	230	295	230	360	520	865
Distance ① (mm)	200	200	200	230	295	490	390	520	615

Q30.195.086.WD

Shelves/number

	4	3
Distance ④ (mm)	370	
Distance ③ (mm)	285	400
Distance ② (mm)	315	380
Distance ① (mm)	315	410

Standard storage heights of the drawers

	9	8	7	6	5	4+1	4	3	2
--	---	---	---	---	---	-----	---	---	---

Standard storage heights of the shelves

	4	3
Distance ④ (mm)	370	
Distance ③ (mm)	285	400
Distance ② (mm)	315	380
Distance ① (mm)	315	410


Q30.195.086.WD

Quantity

Quantity	Type of container	Container size
49	Aerosol cans	400 ml
23	Cans	1 litre
23	Laboratory glass bottles	1 litre
23	Laboratory PE-bottles	1 litre
9	Laboratory glass bottles	2.5 litres

Storage capacities per drawer/type of container

Type of container	Container size
Aerosol cans	400 ml
Cans	1 litre
Laboratory glass bottles	1 litre
Laboratory PE-bottles	1 litre
Laboratory glass bottles	2.5 litres

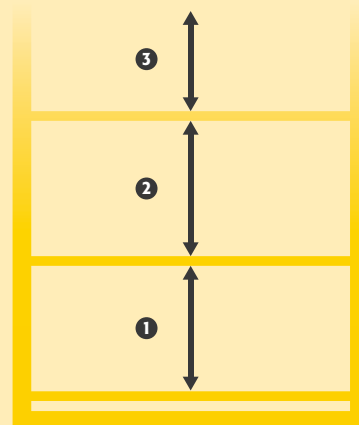



Safe and approved storage of flammable hazardous substances in working areas

Function / construction:

- **robust construction and longevity:** triple hinged door, cabinet body made of melamine resin and flame retardant special plates, safety elements assembled outside the storage compartment for increased protection against corrosion, doors with metal front as well as scratch-and impact-resistant surface, easy to clean
- **easy handling with comfort:** door stops in any position, self-closing in the event of fire
- **no unauthorised use:** door lockable with profile cylinder (integration in an existing locking system possible)
- **easy transport:** integrated transport base with removable, optional base cover or adjusting aids with base cover

- **ventilation:** integrated air ducts ready for connection (DN 75) to a technical exhaust system

Available equipment

- shelves with load capacity of 75 kg
- tray shelves with load capacity of 150 kg
- perforated insert useable as storage area
- bottom collecting sump
- drawers – high storage capacity, optimal use of the available storage capacities
- drawers heavy-duty with load capacity of 100 kg
- metal-free tray shelves and collecting sumps
- storage box
- sump inliner (PE or PP)

Safety storage cabinet

Q-CLASSIC-30 model Q30.195.056

Body colour Anthracite grey, with wing door in Warning yellow (RAL 1004), interior equipment with 3 x shelf, 1 x perforated insert, 1 x bottom collecting sump (sheet steel powder-coated)

Order No. 31002-060-33600

(base cover and extraction unit optional)



Drawers

High storage capacity
Ideal for the storage of small containers



Shelf

Height adjustable (32 mm increments), with 30 mm upstand



Extraction unit

Safe ventilation of safety storage cabinets

plug-in – fast and easy to install

Further information and products from page 397

Order No. 14218



Metal-free interior equipment

Ideal for the storage of flammable, aggressive hazardous substances. Tray shelves with chemical resistant surface lamination and removeable plastic sump, bottom collecting plastic sump.



Recirculating air filter system

Avoid complex wall break-throughs and expensive extraction piping

Further information from page 393


Order No. 25099

**Safety storage cabinet****Q-CLASSIC-30 model Q30.195.056**

Body colour Anthracite grey, with wing door in Traffic red (RAL 3020), interior equipment with 4 x drawer (sheet steel powder-coated)

Order No. 31002-066-33608

(base cover optional)

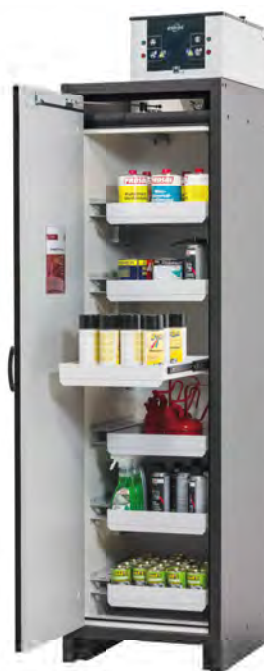**Safety storage cabinet****Q-CLASSIC-30 model Q30.195.056**

Body colour Anthracite grey, with wing door in Light grey (RAL 7035), interior equipment with 6 x drawer (sheet steel powder-coated)

Order No. 31002-061-33610

(recirculating air filter system optional)

7 door colours are available **at no extra cost:**



Order example - how to order correctly



	Model	Information / equipment	Order No.
	Q30.195.056	without interior equipment, left hinged	31002
	Q30.195.056.R	without interior equipment, right hinged	31007

	Body colour	Door colour	Order No.
	Anthracite grey	Warning yellow RAL 1004	060
	Anthracite grey	Light grey RAL 7035	061
	Anthracite grey	Traffic red RAL 3020	066
	Anthracite grey	Gentian blue RAL 5010	067
	Anthracite grey	Reseda green RAL 6011	068
	Anthracite grey	asecos silver	069
	Anthracite grey	Pure white RAL 9010	070

	Interior equipment packages (can only be ordered in combination with the cabinet)	Material	Order No.
			old code new code
	3x shelf, 1x perforated insert, 1x bottom collecting sump (V=22.0L)	sheet steel powder-coated RAL 7035	30076 33600
	4x shelf, 1x perforated insert, 1x bottom collecting sump (V=22.0L)	sheet steel powder-coated RAL 7035	30077 33602
	3x shelf, 1x bottom collecting sump (V=22.0L)	sheet steel powder-coated RAL 7035	30078 33603
	4x shelf, 1x bottom collecting sump (V=22.0L)	sheet steel powder-coated RAL 7035	30079 33604
	3x tray shelf (V=12.0L), 1x perforated insert, 1x bottom collecting sump (V=33.0L)	sheet steel powder-coated RAL 7035	30080 33606
	3x shelf, 1x perforated insert, 1x bottom collecting sump (V=22.0L)	stainless steel	30081 30081
	4x shelf, 1x perforated insert, 1x bottom collecting sump (V=22.0L)	stainless steel	30082 30082
	3x shelf, 1x bottom collecting sump (V=22.0L)	stainless steel	30083 30083
	4x shelf, 1x bottom collecting sump (V=22.0L)	stainless steel	30084 30084
	3x tray shelf (V=12.0L), 1x perforated insert, 1x bottom collecting sump (V=33.0L)	stainless steel	30085 30085
	3x tray shelf (V=11.0L), 1x bottom collecting sump (V=11.0L)	PP grey	HIGHLIGHT 30086 30086
	4x tray shelf (V=11.0L), 1x bottom collecting sump (V=11.0L)	PP grey	HIGHLIGHT 30087 30087
	4x drawer (V=5.0L)	sheet steel powder-coated RAL 7035	30088 33608
	5x drawer (V=5.0L)	sheet steel powder-coated RAL 7035	30089 33609
	6x drawer (V=5.0L)	sheet steel powder-coated RAL 7035	30090 33610
	3x drawer (V=5.0L), 1x drawer (disposal) (V=30.0L)	sheet steel powder-coated RAL 7035	30091 33612
	4x drawer (V=5.0L), 1x drawer (disposal) (V=30.0L)	sheet steel powder-coated RAL 7035	30092 33613
	4x drawer (V=5.0L)	stainless steel	30093 33616
	5x drawer (V=5.0L)	stainless steel	30094 33617
	6x drawer (V=5.0L)	stainless steel	30095 33618
	3x drawer (V=5.0L), 1x drawer (disposal) (V=30.0L)	stainless steel	30096 33619
	4x drawer (V=5.0L), 1x drawer (disposal) (V=30.0L)	stainless steel	30097 33620



Individual configuration of model Q30.195.056



Model	Information / equipment	W x D x H (mm)	Order No.
Q30.195.056	without interior equipment, left hinged	564 x 620 x 1947	31002
Q30.195.056.R	without interior equipment, right hinged	564 x 620 x 1947	31007



Body colour	Door colour	Order No.
Anthracite grey	Warning yellow RAL 1004	060 4
Anthracite grey	Light grey RAL 7035	061 4
Anthracite grey	Traffic red RAL 3020	066 4
Anthracite grey	Gentian blue RAL 5010	067 4
Anthracite grey	Reseda green RAL 6011	068 4
Anthracite grey	asecos silver	069 4
Anthracite grey	Pure white RAL 9010	070 4



Interior equipment	Material	W x D x H (mm)	 ¹⁾	 ²⁾	 ³⁾	Order No.	
						old code	new code
 shelf	sheet steel powder-coated RAL 7035	440 x 432 x 30*	75.00			30098	30098 2
	stainless steel	440 x 432 x 30*	75.00			30102	30102 2
 tray shelf	sheet steel powder-coated RAL 7035	427 x 442 x 70*	150.00	12.00	10.91	30099	30099 2
	stainless steel	427 x 442 x 70*	150.00	12.00	10.91	30103	30103 6
 tray shelf, metal-free	PP grey	424 x 472 x 67*	50.00	13.00	11.82	30106	30106 2
 perforated insert	sheet steel powder-coated RAL 7035	420 x 439 x 60	75.00			5555	5555 2
	stainless steel	425 x 439 x 60	75.00			5776	5776 2
 bottom collecting sump	sheet steel powder-coated RAL 7035	442 x 507 x 127		22.00	20.00	30100	33594 2
	stainless steel	418 x 433 x 120*		22.00	20.00	30104	30104 2
 bottom collecting sump, metal-free	PP grey	422 x 472 x 67*		12.50	11.36	25424	25424 2
 bottom collecting sump, height = 180 mm	sheet steel powder-coated RAL 7035	409 x 474 x 171*		33.00	30.00	30101	33605 2
 bottom collecting sump, height = 170 mm	stainless steel	439 x 454 x 170*		33.00	30.00	30105	30105 6
 drawer	sheet steel powder-coated RAL 7035	345 x 500 x 50*	25.00	5.00	4.55	30108	33607 2
	stainless steel	345 x 500 x 50*	25.00	5.00	4.55	30110	33614 2
 drawer (disposal)	sheet steel powder-coated RAL 7035	329 x 462 x 200*	50.00	30.00		30109	30109 2
	stainless steel	329 x 462 x 150*	50.00	30.00		30111	30111 2
 storage box, left hinged	sheet steel powder-coated RAL 7035	380 x 420 x 360				30029	30029 2
 storage box, right hinged	sheet steel powder-coated RAL 7035	380 x 420 x 360				30030	30030 2
 sump inliner, for drawer disposal	PP grey	318 x 449 x 66*				18472	18472 2
 sump inliner	PP grey	353 x 406 x 88*				9706	9706 2
 sump inliner, for drawer disposal	PP black	313 x 413 x 171*				18474	18474 2
 anti-slip rubber mat, for drawer disposal	black	495 x 340 x 1				8719	8719 2
 pipe lead-through	sheet steel powder-coated RAL 7035	182 x 104 x 48				30035	30035 2
 base cover, with adjusting aids	panel melamine resin-coated RAL 7016	488 x NaN				31004	31004 2
 base cover	panel melamine resin-coated RAL 7016					31923	31923 2
 Q-mover, 2 pieces	steel powder-coated textured RAL 7016	600 x 140 x 130	450.00			30037	30037 2
 box for chemicals	PP blue	410 x 300 x 260		18.00	16.36	30038	30038 2




Accessories	Order No.
recirculating air filter system with exhaust air monitoring	25099 2
extraction unit without exhaust air monitoring	14218 2

* usable storage space W x D (mm)


1) load capacity [kg], 2) retention capacity [ltr], 3) maximum container volume in the cabinet [ltr]




Front view

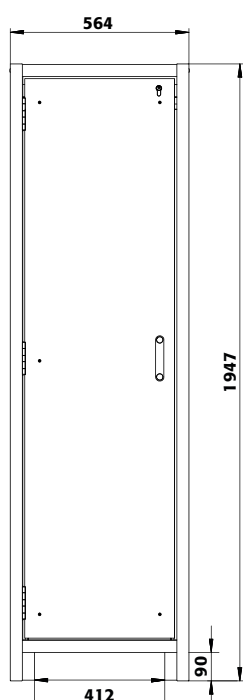

Side view

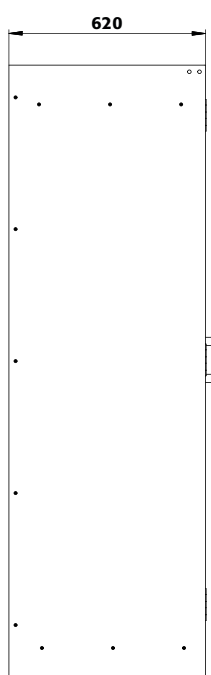

Selectional view

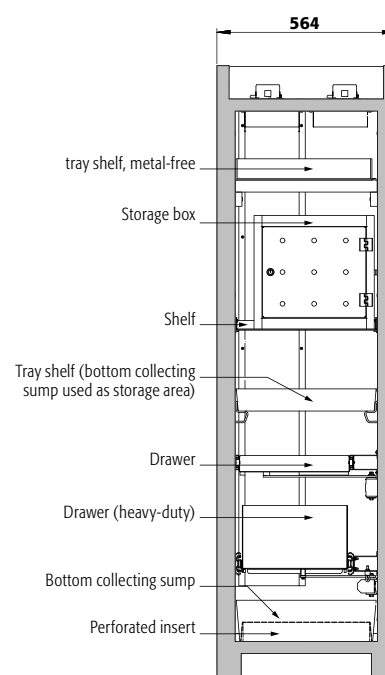

Technical data


Q30.195.056

External dimensions W x D x H	mm	564 x 620 x 1947
Internal dimensions W x D x H	mm	450 x 519 x 1626
Weight without interior equip.	kg	146
Air change 10 times	m ³ /h	4
Diff. pressure (air change 10 times)	Pa	10
Diff. pressure (air change 10 times EK5/AK4)	Pa	10
Maximum load	kg	600
Distributed load	kg/m ²	816

Transport base

Entry height transport base	mm	90
Entry width transport base	mm	412

Top view

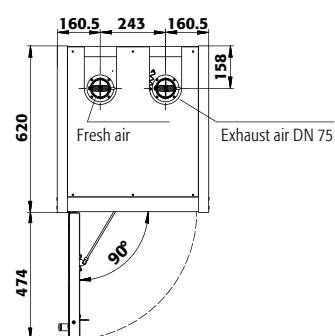

Q30.195.056

Drawers/number

	9	8	7	6	5	4	3	2
Distance ⑨ (mm)	125							
Distance ⑧ (mm)	135	125						
Distance ⑦ (mm)	135	170	155					
Distance ⑥ (mm)	135	170	200	220				
Distance ⑤ (mm)	135	170	200	230	220			
Distance ④ (mm)	135	170	200	230	295	315		
Distance ③ (mm)	135	170	200	230	295	360	410	
Distance ② (mm)	170	170	200	230	295	360	520	860
Distance ① (mm)	200	200	200	230	295	390	520	615

Q30.195.056


Shelves/number

	4	3
Distance ④ (mm)	370	
Distance ③ (mm)	285	400
Distance ② (mm)	315	380
Distance ① (mm)	315	410

Standard storage heights of the shelves

Storage capacities per drawer/type of container

Quantity	Type of container	Container size
30	Aerosol cans	400 ml
12	Cans	1 litre
12	Laboratory glass bottles	1 litre
12	Laboratory PE-bottles	1 litre
6	Laboratory glass bottles	2.5 litres

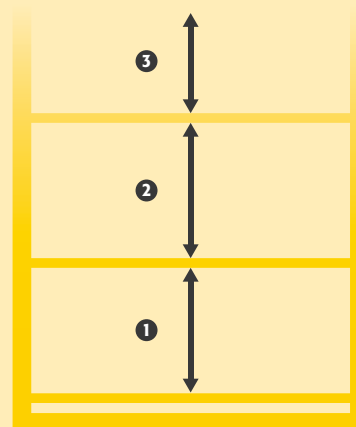



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.