

# BP BULK POWDER ENCLOSURE



**CTS Bulk Powder enclosure systems are designed to provide safe weighing and manipulation environments for powder transfer and dispensing operations.**

Constructed in stainless steel, BP enclosure systems are very robust and easy to clean. The unique airflow design around the drum sampling area provides a high level of containment for the handling of kilogram quantities of active drug compounds. CTS Bulk Powder systems provide an effective alternative to a built in downflow booth system.


- Ergonomic design - Easy to use
- Dual sensor airflow alarm
- Mobile bench mounted solution
- Safe change HEPA filtration
- No facilities extract required
- Custom built sizes available

# BP BULK POWDER ENCLOSURE

## SYSTEM SPECIFICATION



## CONTAMINATION CONTROL

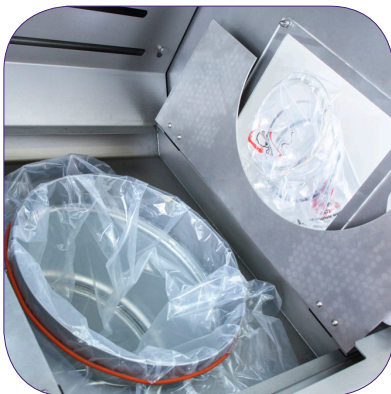

## AIR FLOW & FILTRATION



## CONSTRUCTION

- 316 Stainless steel main body
- 0.5ra Pharmaceutical grade internal surface finish
- 304 Stainless mobile frame
- 8mm Optically clear cast acrylic hinged access door

## DIMENSIONS (mm) - (WxDxH)

External:	Internal:
1200 x 1000 x 1650	1180 x 640 x 720
1500 x 1000 x 1650	1480 x 640 x 720
1800 x 1000 x 1650	1780 x 640 x 720

- Custom built solutions available.

## POWER REQUIREMENT (SINGLE PHASE)

- 230v/50Hz/13A Service (EU)
- 110v/60Hz/20A Service (USA)

## FACILITIES REQUIREMENT

- Re-circulatory design, can be stand-alone operation or connected direct to house extract

## CONTAMINATION CONTROL

Additional extract is present around the drum access area providing additional protection from airborne contaminants released when opening drum bags.

When the sample drum is docked, the user is presented with a level and easily cleanable surface.

## FILTRATION

- Safe-change HEPA filter
- Modular air handling filtration unit
- Min 99.995% efficiency at 0.3 microns (H14)
- Activated Carbon filtration available

## AIRFLOW MONITORING

- Continual airflow monitoring using a dual sensor hot wire anemometer complete with Audible / Visual alarm



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.