

UV-visible spectrophotometer

**SP-V1000/SP-UV1000/
SP-V1100/ SP-UV1100**



Features

- The shape is professional, concise, friendly and mellow, and the appearance design conforms to ergonomics. The use of high-quality materials, fine craftsmanship, soft texture, showing a good user experience
- Attention to human-device experience, friendly interactive interface, easy for users to operate and detect
- Up to 128 sets of standard curves and 4000 test data can be saved
- Test data is saved in real time to avoid data loss in case of sudden power failure.
- Automatic wavelength calibration and dark current correction
- It can automatically switch the light source and turn off the corresponding light source, so as to prolong the service life of the light source to the greatest extent.
- The instrument has rich functions and can realize various functional operations, and the host can independently complete photometric measurement and quantitative test, kinetics measurement and other functions
- Can be used with automatic cell holder, film holder, tube rack, constant temperature cell holder, reflectivity accessory, variable optical path sample holder, variable angle solid sample holder and other special accessories
- The computer software is designed with single interface and multiple functions, which has friendly interface, simple operation, rich and powerful functions, and fast response speed. Real-time monitoring data change. The unique shortcut toolbar makes the experimental process clear at a glance, which greatly facilitates the use. And can be extended more application functions, such as wavelength scanning, time scanning (kinetics), multi-wavelength scanning, and DNA/protein detection analysis.


Applications

The spectrophotometer produced by our company refers to the absorbance or transmittance of a substance at a specific wavelength or within a certain wavelength range, carry out qualitative or quantitative analysis on that substance, as well as kinetic experiments, protein or DNA detection, DNA and RNA quantitative analysis, protein quantitative analysis and other applications. It can be widely used in medicine and health, clinical testing, biochemistry, petrochemical industry, environmental protection testing, food hygiene and quality control. It can also be used as a teaching demonstration and experimental instrument for related courses in colleges and universities.


Photometry

546.0nm 0.001A		
No.	WL.	Abs.
1	230.0	0.001
2	340.0	0.000
3	450.0	0.002
4	540.0	0.000
5	620.0	0.003

Quantitative

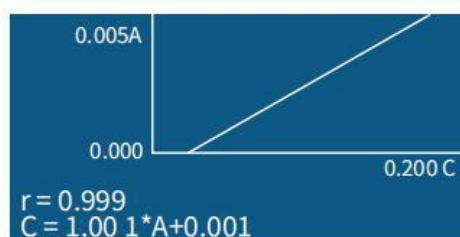

Kinetics

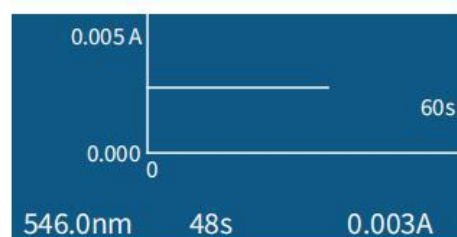



Model	SP-V1000	SP-UV1000	SP-V1100	SP-UV1100
Optical System	Single Beam	Single Beam	Single Beam	Single Beam
Light Source	Tungsten lamp	Tungsten lamp & Deuterium lamp	Tungsten lamp	Tungsten lamp & Deuterium lamp
Spectral Bandwidth	4.0nm	4.0nm	2.0nm	2.0nm
Wavelength Range	325-1000nm	200-1000nm	320-1100nm	190-1100nm
Wavelength Accuracy	±2nm	±2nm	±0.5nm	±0.5nm
Wavelength Repeatability	≤1nm	≤1nm	≤0.3nm	≤0.3nm
Photometric Range	-0.3 - 3.0 A, 0 - 200%T	-0.3 - 3.0 A, 0 - 200%T	-0.3 - 3.0 A, 0 - 200%T	-0.3 - 3.0 A, 0 - 200%T
Photometric Accuracy	±0.003 A @ 0.5 A, ±0.006 A @ 1.0 A, ±0.5%T @ 100%T	±0.003 A @ 0.5 A, ±0.006 A @ 1.0 A, ±0.5%T @ 100%T	±0.002 A @ 0.5 A, ±0.004 A @ 1.0 A, ±0.3%T @ 100%T	±0.002 A @ 0.5 A, ±0.004 A @ 1.0 A, ±0.3%T @ 100%T
Photometric Repeatability	≤0.0015 A @ 0.5 A, ≤0.003 A @ 1.0 A, ≤0.25%T @ 100%T	≤0.0015 A @ 0.5 A, ≤0.003 A @ 1.0 A, ≤0.25%T @ 100%T	≤0.001 A @ 0.5 A, ≤0.002 A @ 1.0 A, ≤0.15%T @ 100%T	≤0.001 A @ 0.5 A, ≤0.002 A @ 1.0 A, ≤0.15%T @ 100%T
Stability	≤0.002 A/h @ 500 nm (after preheat 2 hours)	≤0.002 A/h @ 500 nm (after preheat 2 hours)	≤0.002 A/h @ 500 nm (after preheat 2 hours)	≤0.002 A/h @ 500 nm (after preheat 2 hours)
Stray Light	≤0.2%T @ 360nm	≤0.2%T @ 360nm	≤0.05%T @ 360nm	≤0.05%T @ 220 & 360nm
Sample Compartment	4-position, 10mm pathlength cuvette	4-position, 10mm pathlength cuvette	4-position, 10mm pathlength cuvette	4-position, 10mm pathlength cuvette
Display	LCD	LCD	LCD	LCD
Output	USB Port & Parallel Port (Printer)	USB Port & Parallel Port (Printer)	USB Port & Parallel Port (Printer)	USB Port & Parallel Port (Printer)
Power Requirements	110/220VAC, 50/60Hz, 80W	110/220VAC, 50/60Hz, 120W	110/220VAC, 50/60Hz, 80W	110/220VAC, 50/60Hz, 120W
Dimension(L×W×H)	490×360×210mm	490×360×210mm	490×360×210mm	490×360×210mm
Weight	12kg	14kg	12kg	14kg





WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.