

FICHE TECHNIQUE

REFERENCE :

060863

ClearLine® Pipette Controller

Description

- Designed for cordless work with 0.5-100 ml glass or plastic pipettes.
- Lightweight handle, smooth pushbuttons and switches ensure effortless pipetting even during extensive use.
- Different operational modes may be selected depending on pipetting volume and viscosity of liquid.
- 2 Liquid aspiration speeds HIGH or LOW, while dispensing can be by gravity (GRAV) or supported by the pump (BLOW).
- PTFE 0.2µm filter protects the unit against overfilling.
- The powerful re-chargeable Ni-MH battery allows many hours of continuous work.
- Supplied with a charging stand & wall mount for convenient charging and storage.
- Rest wings allows practical storage during operation.


Technical specification

- Volume Range: 0.5–100 ml
- Working/Recharging Temperature 10–35 °C
- Storage Temperature -25 to 50 °C
- NiMH Battery AA type 1300 mAh, 2*1 2,4 V
- Autoclavable nosepiece at 121 °C

Packaging:

- White glossy box with label or white box with sleeve
- Instruction manual
- Universal charger 9V with set of plugs: EU, US, UK, AU
- Wall hanger
- • Charging stand
- • PTFE 0.2µm filter

Product Attributes

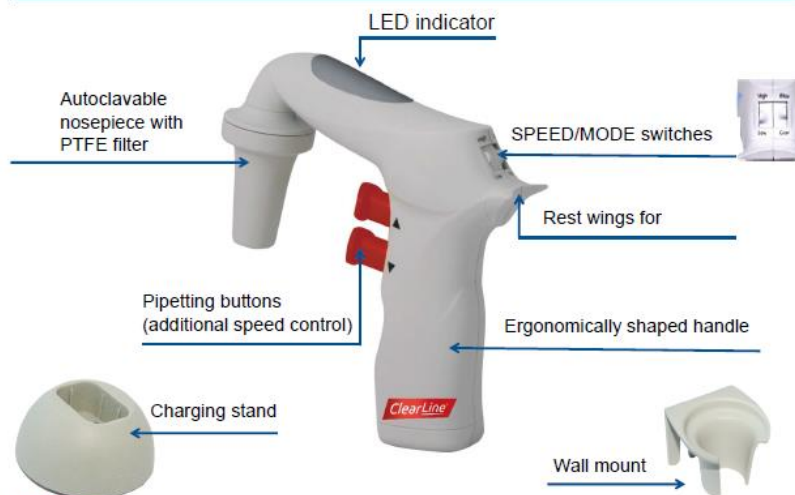



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.