

# Combined Refrigerator and Freezer Unit



HYCD-469A

HYCD-469

## Scope of Application

This product is designed to store frozen plasma, biomaterials, vaccines, electronic devices and low temperature tests of special materials, in blood stations, hospitals, epidemic prevention stations, universities, scientific research institutes, electronic chemical industry and other enterprise laboratories.

## Innovative Design

- ◎ Hydrocarbon Energy Saving
- ◎ Environmentally Friendly
- ◎ Safe, Secure and Reliable
- ◎ Intelligent Control

## Product Advantages



### Frequency Conversion and Energy Saving

Frequency conversion high efficiency compressor, optimized refrigeration system, more energy saving.



### Green Environmental Protection

- Hydrocarbon compressor, hydrocarbon refrigerant, completely green refrigeration system.
- Cyclopentane foamed polyurethane insulation layer, the whole series conform to the European ROHS directive material, green environmental protection.



### Large Space, Good Uniformity, Accurate Temperature Control

- Built-in evaporator, large space in the box, shelf/partition adjustable.
- Temperature uniformity of refrigerator  $\pm 3^{\circ}\text{C}$ ; Control accuracy  $0.1^{\circ}\text{C}$ .
- Freezer temperature uniformity  $\pm 5^{\circ}\text{C}$ ; Control accuracy  $1^{\circ}\text{C}$ .



### Stainless Steel Inner for Easy Cleaning

Stainless steel inner tank is corrosion resistant and easy to clean.



### Microprocessor Control System

- Microprocessor temperature controller;  
Refrigerator adjustable temperature  $2\sim 8^{\circ}\text{C}$ , display accuracy at  $0.1^{\circ}\text{C}$ ;  
Freezer adjustable temperature  $-10\sim -40^{\circ}\text{C}$ , display accuracy at  $1^{\circ}\text{C}$ .
- Multiple alarm functions including high/low temperature alarm, sensor error alarm, power failure alarm, low battery alarm, door ajar alarm and high ambient temperature alarm.
- Two types of alarm indications: audible buzzer and visible flashing light. Alarm lasts for more than 48 hours after power failure.
- Standard USB & optional printer.

## Ergonomic Design



Equipped with padlocks in both refrigerator and freezer to ensure storage safety



- Equipped with electrically heated glass to prevent condensation effectively (HYCD-469A)
- Equipped with overcurrent protection to ensure usage safety



Two test holes for both refrigerator and freezer as standard, allow independent testing of cabinet temperature



Refrigerator and freezer are equipped with self-closing hinges to avoid door ajar (less than  $90^{\circ}$  self-closing)

# Combined Refrigerator and Freezer Unit

## Product Parts Diagram



- Microcomputer controller with dual display and separate control
- Microprocessor temperature control accuracy 0.1°C (refrigerator) and 1°C (freezer)

USB interface for easier data transfer

Dual fan system providing high refrigeration efficiency with a better temperature uniformity and stability

Equipped with LED lights in the refrigerator

Equipped with three adjustable dip plastic shelves in the refrigerator

Glass door (HYCD-469A)

Equipped with padlocks in both refrigerator and freezer to ensure storage safety

Stainless steel tank


Equipped with three adjustable stainless steel shelves in the freezer, with three optional plastic drawers for different needs

Four casters and two levelling legs for ease of maneuvering and positioning




## 32°C Ambient Temperature Curve


### Refrigerator Cool Down Curve

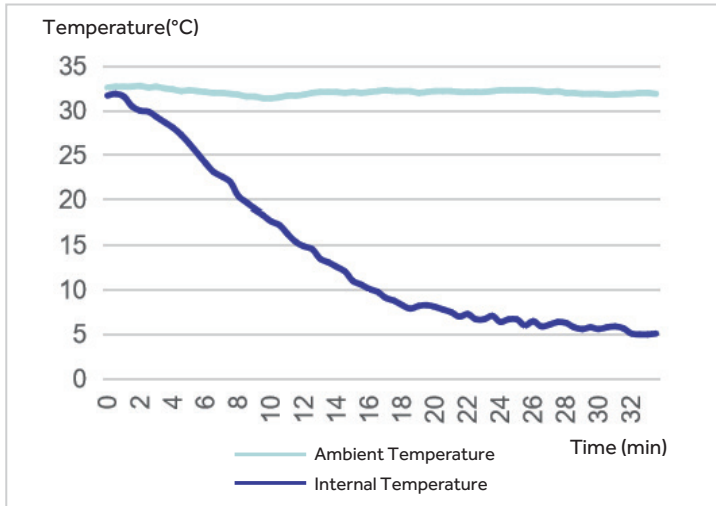


### Refrigerator Warm Up Curve

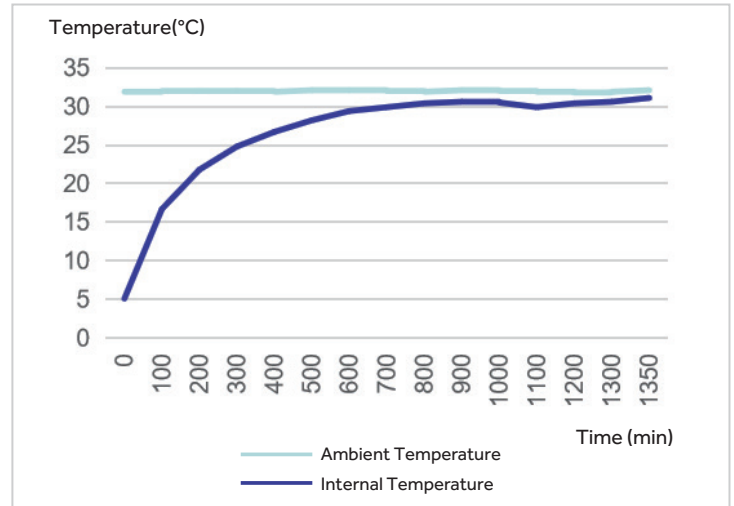


### Freezer Cool Down Curve

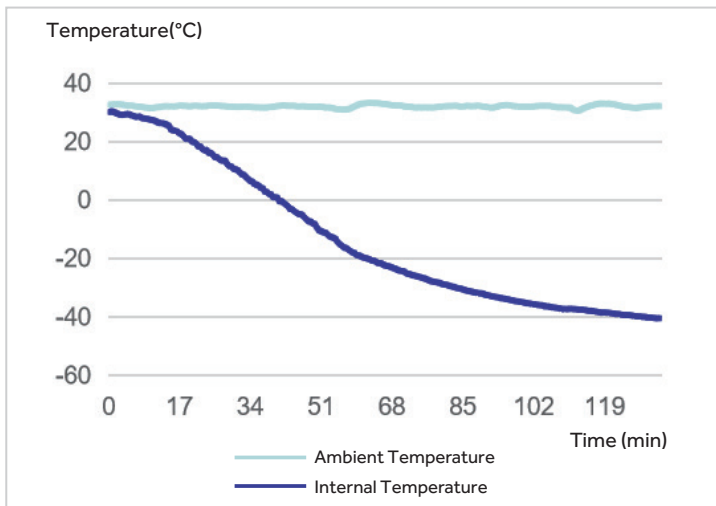

### Freezer Warm Up Curve

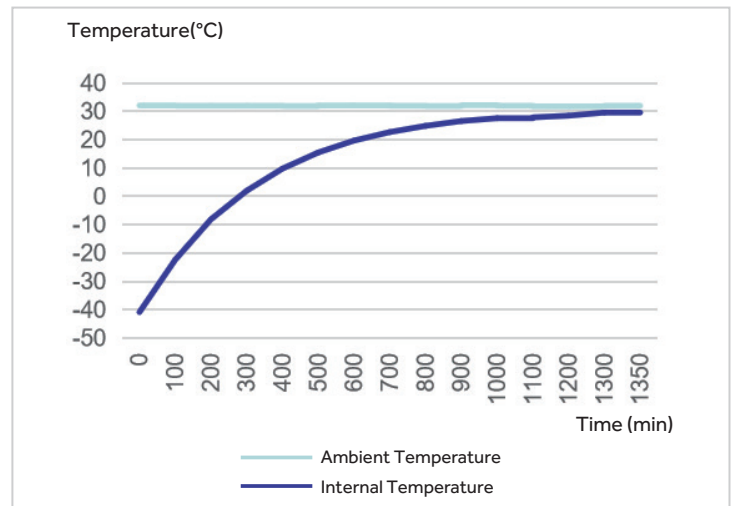

## Dimensions

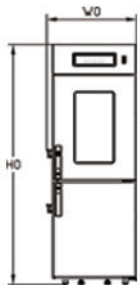

W0-Total Width  
W1-Refrigerator Width  
W2-Freezer Width  
D0-Total Depth  
D1-Refrigerator Depth  
D2-Freezer Depth  
H0-Total Height  
H1-Refrigerator Height  
H2-Freezer Height

| W0  | W1  | W2  | D0  | D1  | D2  | H0   | H1  | H2  |
|-----|-----|-----|-----|-----|-----|------|-----|-----|
| 738 | 593 | 553 | 834 | 585 | 545 | 1975 | 800 | 660 |

## Specifications



| Model           |                               |     | HYCD-469                                                                         | HYCD-469A                                                   |
|-----------------|-------------------------------|-----|----------------------------------------------------------------------------------|-------------------------------------------------------------|
| Technical Data  | Cabinet Type                  |     | Upright                                                                          | Upright                                                     |
|                 | Climate Class                 |     | N (16~32°C)                                                                      | N (16~32°C)                                                 |
|                 | Cooling Type                  |     | Refrigorator: Forced air cooling<br>Freezer: Direct cooling                      | Refrigorator: Forced air cooling<br>Freezer: Direct cooling |
|                 | Defrost Mode                  |     | Refrigorator: Auto / Freezer: Manual                                             | Refrigorator: Auto / Freezer: Manual                        |
|                 | Refrigerant                   |     | Refrigorator: R600a / Freezer: R290                                              | Refrigorator: R600a / Freezer: R290                         |
|                 | Sound Level (dB(A))           |     | 45                                                                               | 45                                                          |
| Performance     | Temperature Range (°C)        |     | Refrigorator: 2~8 / Freezer: -10~-40                                             | Refrigorator: 2~8 / Freezer: -10~-40                        |
| Control         | Controller                    |     | Microprocessor                                                                   | Microprocessor                                              |
|                 | Display                       |     | LED                                                                              | LED                                                         |
| Electrical Data | Power Supply (V/Hz)           |     | 220~240/50/60                                                                    | 220~240/50/60                                               |
|                 | Power (W)                     |     | 640                                                                              | 640                                                         |
|                 | Electrical Current (A)        |     | 3                                                                                | 3                                                           |
| Dimensions      | Capacity (L/Cu.Ft)            |     | Refrigorator: 269/9.5; Freezer: 200/7.1                                          | Refrigorator: 269/9.5; Freezer: 200/7.1                     |
|                 | Net/Gross Weight (approx)     | kg  | 176/207                                                                          | 176/207                                                     |
|                 |                               | lbs | 388.0/456.4                                                                      | 388.0/456.4                                                 |
|                 | Interior Dimensions (W*D*H)   | mm  | Refrigorator: 593*585*800<br>Freezer: 553*545*660                                | Refrigorator: 593*585*800<br>Freezer: 553*545*660           |
|                 |                               | in  | Refrigorator: 23.3*23.0*31.5<br>Freezer: 21.8*21.5*26.0                          | Refrigorator: 23.3*23.0*31.5<br>Freezer: 21.8*21.5*26.0     |
|                 | Exterior Dimensions (W*D*H)   | mm  | 738*834*1975                                                                     | 738*834*1975                                                |
|                 |                               | in  | 29.1*32.8*77.8                                                                   | 29.1*32.8*77.8                                              |
|                 | Packing Dimensions (W*D*H)    | mm  | 800*875*2100                                                                     | 800*875*2100                                                |
|                 |                               | in  | 31.5*34.4*82.7                                                                   | 31.5*34.4*82.7                                              |
|                 | Container Load (20'/40'/40'H) |     | 14/28/28                                                                         | 14/28/28                                                    |
| Alarms          | High/Low Temperature          |     | Y                                                                                | Y                                                           |
|                 | Sensor Error                  |     | Y                                                                                | Y                                                           |
|                 | Power Failure                 |     | Y                                                                                | Y                                                           |
|                 | Door Ajar                     |     | Y                                                                                | Y                                                           |
|                 | Remote Alarm                  |     | Y                                                                                | Y                                                           |
|                 | High Ambient Temperature      |     | Y                                                                                | Y                                                           |
|                 | Low Battery                   |     | Y                                                                                | Y                                                           |
| Accessories     | Caster                        |     | Y                                                                                | Y                                                           |
|                 | Foot                          |     | Y                                                                                | Y                                                           |
|                 | Porthole                      |     | Y                                                                                | Y                                                           |
|                 | Shelves/Drawers               |     | Refrigorator: 3/0<br>Freezer: 3/0 (Max 3 Drawers optional, without sliding rack) |                                                             |
|                 | USB Interface                 |     | Y                                                                                | Y                                                           |
| Others          | Certification                 |     | CE                                                                               | CE                                                          |



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.