



# 90 Gallon, 2 Shelves, 2 Doors, Self Close, Flammable Cabinet, Sure-Grip® EX, Yellow - 899020

## Product Images

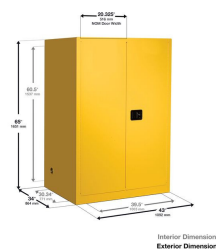

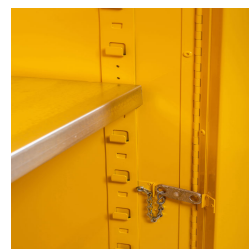

## Short Description

---

- Designed for the safe storage of chemicals and flammable liquids
- Exterior dimensions of 65 in H x 43 in W x 34 in D (1651 mm H x 1092 mm W x 864 mm D)
- Storage capacity of 90 gal (341 L) of flammable liquids
- Fire-resistant, leak-tested design
- Complies with OSHA, NFPA 30, NFPA 1 and IFC
- Fire tested and FM-approved

## Description

---

One of the leading causes of industrial fires is the improper storage and handling of flammable liquids. The Justrite Standard Sure-Grip® EX Flammable Safety Cabinet with Self-Close Doors (model 899020) is designed to safely store flammable liquids. Its exterior dimensions are 65 in H x 43 in W x 34 in D (1651 mm H x 1092 mm W x 864 mm D). These yellow cabinets can store 90 gal (341 L)\* of flammable liquids.

This model is constructed with a sturdy 18-gauge (1-mm) thick steel double-wall, with 1.5 in (38 mm) of insulating air space for fire resistance. The fully welded, not riveted, construction holds the squareness of the unit for a longer life and offers greater protection in a fire since air gaps are reduced.

High-performance, self-closing doors open to 140 degrees. A patented, concealed self-close mechanism provides obstruction-free access to shelf space. The doors self-index to guarantee they will close in tight sequence. Continuous piano hinges ensure the doors close smoothly every time.

When used properly as designed, safety cabinets with self-closing doors shut automatically in case of a fire. The self-closing doors, under tension, are held open by using fusible links. Under fire conditions, when the ambient temperature reaches 165 F (74 C), the links melt and the safety cabinet's self-closing doors automatically swing shut. The self-latching system engages with three stainless steel bullet-slams located at the top, bottom, and center locations of the door. For replacement fusible links, order model 29988.

For added security, each cabinet comes with an innovative U-Loc™ handle that does not require manual rotation to engage the mandatory latch. Cabinets include a cylinder lock and key set, or you can choose to use your own padlock. Trilingual Haz-Alert™ reflective warning labels that read "FLAMMABLE – KEEP FIRE AWAY" are highly visible under fire conditions or during power outages when illuminated with a flashlight.

Model 899020 features two patented, galvanized steel SpillSlope® shelves designed to direct spills to the back and bottom into a 2-in (51-mm) leak-tested sump that meets EPA requirements. Adjustable on 3-in (76-mm) centers, all SpillSlope shelves meet ANSI MH28.1 standards with a 350-lb (159-kg) safe allowable load. Welded shelf hangers interlock with the shelf to offer maximum stability with no flimsy brackets to lose. For replacement shelves, use model 29945.

The cabinet also features dual vents with flame arresters, four adjustable self-leveling feet, and a built-in grounding connector. Its durable lead-free epoxy/polyester powder-coat finish provides excellent chemical resistance inside and out, as well as minimizes the effects of corrosion and humidity.

Justrite supplies more than 300 safety cabinet models ranging from 4- to 115-gal (15- to 435-L), to meet a versatile range of safety, compliance, and storage requirements.

All Justrite cabinets carry a 10-year warranty. Model 899020 is tested and approved by FM Approvals. It complies with OSHA, NFPA 30, NFPA 1, and IFC requirements.

## Safety Considerations for Standard Safety Cabinets

The requirements for handling and storing chemicals are complex and depend on how flammable the liquid is. There is no room for error, so it is critical to choose a safety cabinet that has been independently tested to meet the standards of your industry and local jurisdictions. Not just any safety cabinet will do. Justrite flammable cabinets are best in class and are trusted within a wide range of industrial settings across the world, including laboratories, manufacturing, oil refineries, chemical processing plants, and more.

For more information, please consult our workplace safety article about the [requirements for flammable liquid storage](#).

\*OSHA 29 CFR 1910.106(d)(3) states that not more than 60 gallons (227L) of Category 1, 2, or 3 flammable liquids, nor more than 120 gallons (454L) of Category 4 flammable liquids may be stored in a safety cabinet.

## Specifications

|                                       |                               |
|---------------------------------------|-------------------------------|
| Model No                              | 899020                        |
| Badge                                 | Best Seller                   |
| UPC                                   | 697841114845                  |
| International Model No.               | 8990201                       |
| Color                                 | Yellow                        |
| Pallet Quantity                       | 1                             |
| Material Specifications               | 18-Gauge Steel                |
| FM Approval                           | Yes                           |
| OSHA Compliance                       | Yes                           |
| Complies with International Fire Code | Yes                           |
| Complies with NFPA Code 30            | Yes                           |
| Complies with NFPA 1 Fire-Code        | Yes                           |
| Capacity                              | 90 gallon                     |
| Liter Capacity                        | 341                           |
| Number of Doors                       | 2                             |
| Door Type                             | Self-close                    |
| Number of Shelves                     | 2                             |
| Adjustable Shelves                    | Yes                           |
| Dimensions, Exterior                  | 65" H x 43" W x 34" D         |
|                                       | 1651mm H x 1092mm W x 864mm D |
| Dimensions, Interior                  | 60.5" H x 39.5" W x 30.3" D   |
|                                       | 1537mm H x 1003mm W x 771mm D |
| Height (Inches)                       | 65                            |
| Width (Inches)                        | 43                            |
| Depth                                 | 34                            |
| Net Weight, lbs                       | 381.0                         |

Net Weight, kg

172.8



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.