

Pharmacy Plus Bluetooth Enabled Refrigerator

120L Pharmacy Plus Undercounter Fridge with Glass Door

Product code PPGR120BT

SKU code 444411884




Key features

- ✓ Temperature range of +2°C to +8°C
- ✓ Audible and visual temperature alarms
- ✓ Solid and glass door options
- ✓ Integrated lock, supplied with 2 keys for added security



Bluetooth Connected Controller

Our controller is Bluetooth enabled and allows for simple, contactless monitoring of temperature and alarm history, via the official Lec Medical app. Battery back-up, provides continuous temperature monitoring in the event of power outage.

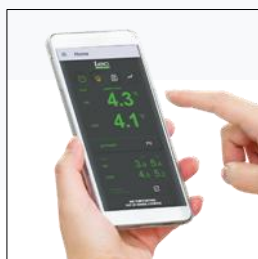

Lec Medical Smartphone App

Contactless local monitoring of air and load temperatures via smart phone app. Downloadable temperature report and alarm history.



Intelligent Fan Management

Smart fan optimisation aides in temperature stability following door openings, minimising fan activity when door is open to reduce the rate of temperature increase.



Dual Probe Temperature Monitoring

Our dual probe system offers greater temperature accuracy within your appliance, simultaneously measuring both air and product temperatures.



Flexible Internal Storage

Multiple shelf locations provides flexible storage solutions.

Specification: 120L Pharmacy Plus Undercounter Fridge with Glass Door

Product Information		Door Specification	
Brand	Lec Medical	Door Finish	PVC / Glass
Category	Pharmacy	Door Type	Glass
Sub-Brand	Pharmacy Plus	Door Style	Hinged reversible
SKU Code	444411884	Door Liner	n/a
EAN Code	5052263118843	Number of Door Trays	n/a
GDMA Code	PPGRI20BT	Number of Door Tray Positions	n/a
Description	Under Counter	Door Open Depth (mm)	1035mm
		Replacement Door Seal	Yes
Product Features		Technical Information	
Gross Capacity (litres)	120L	Noise Level (dB)	43.0dB
Net Capacity (litres)	112L	Colour	White/green
Temperature Range	+2°C to +8°C	External Finish	Painted steel
Open Door Sensor	Yes	Internal Finish	HIPS
Key Lock	Yes, 2 Keys	Type of plastic	SRL800 (HIPS/ABS mix)
Digital Lock	No	Ambient Temperature Range	16°C to 32°C
Porthole	Yes	Refrigerant Type / Weight (g)	R600a/26g
Porthole Size (Ø Diameter)	16mm	CFC Free	Yes
Porthole Location	Rear Bottom (Right)	HFC Free	Yes
Controller Type	Digital BT, External	HCFC Free	Yes
Temperature Display on Panel	Yes	Ammonia Free	Yes
Temperature Alarm	Yes	Insulation Thickness (mm)	40mm
Remote Alarm Contacts	Yes	Insulation Type	High Density Polyurethane Foam
Audible Alarms	Yes	Blowing Agent	Cyclopentane
Visual Alarms	Yes	Certification	CE/UKCA
Interior Light	Yes	ATEX Zone Classification	n/a
Light Type	LED		
Light Activation	Door Switch / App		
Fan - Type	Yes - Axial		
Spark Free	No		
Defrost setting	Automatic		
Battery Back Up/Hours	Yes/48hrs		
Bluetooth	Yes		
NFC	Yes		
WIFI	No		
GSM-4G	No		
Number of Probes	2		
Probe Type	Air & Load		
Anti Microbial Handle	Yes		
Internal Specification		Power/Installation	
Number of Shelves	3	UK 3 Pin Plug Power Cord	Yes
Number of Shelf Positions	7	EU Schuko 2 pin plug Power Cord	Yes
Shelf Size 1 (mm) (WxDxH)	400mm x 390mm x 6mm	AU Type I Plug Power Cord	No
Shelf Size 2 (mm) (WxDxH)	400mm x 285mm x 6mm	Voltage Supply (V)	220-240V
Shelf Material	Polythene Coated Wire	Frequency (Hz)	50/60Hz
Number of Drawers	n/a	Max Power (W)	86.0W
Number of Drawer Positions	n/a	Max Current Rating (A)	0.60A
Drawer Size 1 (mm) (WxDxH)	n/a	Energy Consumption (kWh/24h)	0.52kWh/24h
Drawer Size 2 (mm) (WxDxH)	n/a	Energy Consumption AE (kWh)	199kWh
Drawer Material	n/a	Heat Output Estimated at 22°C (kJ)	1049kJ/h
Number of Baskets	1	Adjustable Feet (Number & position)	2 adjustable (front)
Basket Material	Polythene Coated Wire		
Basket Size (mm) (WxDxH)	375mm x 295mm x195mm		
		Weights & Dimensions	
		Unit External Width (mm)	475mm
		Unit External Depth (mm)	595mm
		Unit External Height (mm)	845mm
		Unit External Depth Inc Handle (mm)	635mm
		Internal Dimensions (mm) (WxDxH)	385mm x 495mm x 700mm
		Product Weight (kg)	34.78Kg
		Packed Width (mm)	505mm
		Packed Depth (mm)	660mm
		Packed Height (mm)	865mm
		Packed Weight (kg)	38.60Kg
		Country of Origin	China



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.