


RLDG0119 Pharmacy & Vaccine Refrigerator

Part Number	RLDG0119
Capacity (Litres)	36
Style	Bench Top
Exterior Dimensions (HxWxD)	540 x 450 x 510mm
Interior Dimensions (HxWxD)	392 x 331 x 226mm
Operating Temperature	+5°C
Temperature Range	+2°C to +8°C
Power Consumption Kwh/24hrs	1.53
Interior Construction	Moulded Plastic
Exterior Construction	Powder Coated Steel
Brackets for Wall Mounting	✓
Shelves	2
Net Weight (kg)	31
Door Style	Glass
Door Hinging	Right (can be reversed at time of manufacture)
Self Closing Door	✓
Key Lock	✓
Digital Temperature Display	✓
Internal Light	✓
High/Low Temperature Alarm	✓
Minimum/Maximum Temperature Recording	✓
Automatic Temperature & Fridge Operation Recording	✓
Door Open Alarm	✓
Mains Failure Alarm	✓
Alarm Battery Back-Up	✓
Second Temperature Probe	✓
Off Cycle Auto Defrost	✓
Fan Assisted Cooling	✓
Remote Alarm Contacts	✓
Access Port	✓
DH Green Book Guidelines Compliant	✓



Powered by the NEW IntelliCold® touch screen controller, this pharmacy refrigerator is packed with extra functions designed to store temperature sensitive pharmaceuticals in optimum conditions. It also features a 1 metre secondary (wandering) temperature probe that can be placed anywhere inside the fridge.

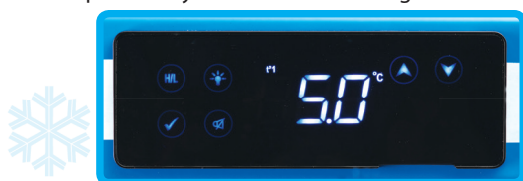

The touch screen controller is completely flat so can be thoroughly cleaned and disinfected as there are no raised buttons to collect dirt and germs. The controller displays, alarms and records minimum and maximum temperatures for both temperature probes. Like a mobile phone touch screen it locks after 30 seconds to ensure that buttons can't be activated accidentally.

For extra security, there is a digital lock version RLDG0119DIGLOCK which solves the problem of generic keys, but whatever the preferred model, for a small additional cost we can provide a manufacturers' calibration certificate from our in house UKAS temperature calibration laboratory No:8898.

We are confident that this pharmacy refrigerator will provide you with many years of fault free service. However, for extra peace of mind your RLDG0119 comes complete with a FREE UK 5 year parts and labour warranty.



While you're registering your warranty, why not purchase a micro SD card and adapter for transferring data recorded by the controller from the fridge to your computer at half price?



The IntelliCold® controller records alarm states, fridge operation, door openings and temperature data for both probes at minute intervals. You can only view and store this data on a computer by purchasing a specially encoded micro SD card from Labcold. This data helps you with record keeping and regulatory inspections, plus if you have an issue, you can establish the exact duration and severity of the event.



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.