



- Extensive Range
- 90+ Years of History
- Extended Aftercare Support
- UKAS Laboratory

### Intellcold® Pharmacy Fridges

## Bench Top RLDG0124DIGLOCK

A next generation pharmacy refrigerator featuring the latest IntelliCold® 24-Range technology, designed for the safe storage of cold chain medicines in line with Public Health England guidance.

An enhanced IntelliCold® touch screen controller delivers accurate timekeeping, intelligent event logging and Bluetooth connectivity for fast access to temperature data via the Labcold Connect Android App. The flat, easy-clean controller locks automatically to prevent accidental changes, while a factory-fitted digital door lock eliminates the risk of lost or shared keys and puts you in control of authorised access. The user code can be changed at any time, with an optional master key available for added reassurance.

Inside, three independent temperature probes provide advanced protection, including a low-temperature safeguard that prevents freezing. A pre-fitted simulated load has been carefully selected to represent the thermal behaviour of commonly stored pharmacy products. This ensures temperature readings reflect true product conditions, providing meaningful alarms when medicines are genuinely at risk; rather than masking excursions with oversized or unrepresentative loads.

Each unit is supplied with a FREE UK 5-year warranty for complete peace of mind. A UKAS temperature calibration certificate (Lab No. 8898) is available at an extra cost if required.

### Specifications

Product Code	RLDG0124DIGLOCK
Gross Capacity	36 litres
Style	Bench Top
UK Warranty	5 Years
External Dimensions (mm) HxWxD	540 x 450 x 510
Internal Dimensions (mm) HxWxD	392 x 450 x 226
Temperature Range	2°C to 8°C
Fridge Set Point	5°C




## Specifications (cont.)

Internal Material	Moulded Plastic
External Material	Powder Coated Steel
Shelves	3
Net Weight	31 Kg
Door Style	Glass
Lock Type	Digital Lock
Digital Temperature Display	Yes
Digital Temperature Control	Yes
High/Low Temperature Alarm	Yes
Min/Max Temperature Recording	Yes
Integral Temp Data Logging	Yes
Automatic Temperature & Fridge Operation Data	Yes
Number of Temperature Probes	3
Internal Light	Yes
Self Closing Door	Yes
Door Open Alarm	Yes
Mains Failure Alarm	Yes
Alarm Battery Backup	Yes
Auto Defrost	Yes
Fan Assisted Cooling	Yes
Remote Alarm Contacts	Yes
Access Porthole	Yes
Access Porthole Size	12mm
Power Consumption Kwh/24hrs	
Ambient Operating Temperature	16°C to 32°C

## Features


- Fitted with the IntelliCold® touch screen controller
- 24/7 automatic temperature logging which is accessible via Bluetooth App
- Fan assisted cooling
- T3 Sensor & Low-Temperature Cut
- Controller and alarm battery back up
- Minimum/maximum temperature recording and display
- Auto defrost
- Self-closing door
- Remote alarm contacts
- FREE five-year labour warranty
- Event based logging
- Digital lock
- Realistic pre-fitted simulated Load
- High/Low temperature alarm
- Door open and power failure alarm
- Access port for additional monitoring equipment
- Internal light
- Brackets for wall mounting

## Models

	Model	Gross Capacity	Door Style	Lock Type
	Bench Top RLDF0119 Diglock	36 litres	Solid	Digital Lock
	Bench Top RLDG0219 Diglock	66 litres	Glass	Digital Lock
	Bench Top RLDF0224	66 Litres	Solid	Key Lock
	Bench Top RLDF0224DIGLOCK	66 Litres	Solid	Digital Lock
	Bench Top RLDF0124	36 litres	Solid	Key lock
	Bench Top RLDG0124	36 litres	Glass	Key Lock
	Bench Top RLDG0224	66 Litres	Glass	Key Lock
	Bench Top RLDG0124DIGLOCK	36 litres	Glass	Digital Lock
	Free Standing RLDF1824	543 Liters	Solid	Key Lock
	Free Standing RLDG1024DIGLOCK	340 Litres	Glass	Digital Lock
	Free Standing RLDG1824	543 Litres	Glass	Key Lock

## Accessories

---

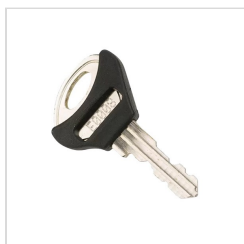

Master Key  
LABLOCK3KEYS



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.