

Pharmacy Plus Bluetooth Enabled Refrigerator

158L Pharmacy Plus Undercounter Fridge with Solid Door


Product code PPSR158BT-UK

SKU code 444411308



Key features

- ✓ Temperature range of +2°C to +8°C
- ✓ Audible and visual temperature alarms
- ✓ Solid and glass door options
- ✓ Integrated lock, supplied with 2 keys for added security

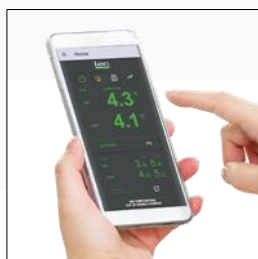

Bluetooth Connected Controller

Our controller is Bluetooth enabled and allows for simple, contactless monitoring of temperature and alarm history, via the official Lec Medical app. Battery back-up, provides continuous temperature monitoring in the event of power outage.

Lec Medical Smartphone App

Contactless local monitoring of air and load temperatures via smart phone app. Downloadable temperature report and alarm history.

Intelligent Fan Management

Smart fan optimisation aides in temperature stability following door openings, minimising fan activity when door is open to reduce the rate of temperature increase.

Dual Probe Temperature Monitoring

Our dual probe system offers greater temperature accuracy within your appliance, simultaneously measuring both air and product temperatures.

Flexible Internal Storage

Multiple shelf locations provides flexible storage solutions.

Product Information	
Brand	Lec Medical
Category	Pharmacy
Sub-Brand	Pharmacy Plus
SKU Code	444411308
EAN Code	5052263113084
GDMA Code	PPSR158BT
Description	Under Counter

Door Specification	
Door Finish	Painted Steel / HIPS*
Door Type	Solid
Door Style	Hinged, Reversible
Door Liner	Flat
Number of Door Trays	n/a
Number of Door Tray Positions	n/a
Door Open Depth (mm)	1150mm
Replacement Door Seal	Yes

Product Features	
Gross Capacity (litres)	158L
Net Capacity (litres)	154L
Temperature Range	+2°C to +8°C
Open Door Sensor	Yes
Key Lock	Yes, 2 Keys
Digital Lock	No
Porthole	Yes
Porthole Size (Ø Diameter)	16mm
Porthole Location	Rear Bottom (Right)
Controller Type	Digital BT, External
Temperature Display on Panel	Yes
Temperature Alarm	Yes
Remote Alarm Contacts	Yes
Audible Alarms	Yes
Visual Alarms	Yes
Interior Light	Yes
Light Type	LED
Light Activation	Door Switch / App
Fan - Type	Yes - Tangential
Spark Free	No
Defrost setting	Auto (Off Cycle)
Battery Back Up/Hours	Yes/48hrs
Bluetooth	Yes
NFC	Yes
WIFI	No
GSM-4G	No
Number of Probes	2
Probe Type	Air & Load
Anti Microbial Handle	Yes

Technical Information	
Noise Level (dB)	401dB
Colour	White/green
External Finish	Painted Steel
Internal Finish	HIPS*
Type of plastic	SRL800 (HIPS/ABS mix)
Ambient Temperature Range	16°C to 32°C
Refrigerant Type / Weight (g)	R600a/28g
CFC Free	Yes
HFC Free	Yes
HCFC Free	Yes
Ammonia Free	Yes
Insulation Thickness (mm)	45mm
Insulation Type	High Density Polyurethane Foam
Blowing Agent	Cyclopentane
Certification	CE/UKCA
ATEX Zone Classification	n/a

Power/Installation	
UK 3 Pin Plug Power Cord	Yes
EU Schuko 2 pin plug Power Cord	Yes
AU Type I Plug Power Cord	No
Voltage Supply (V)	220-240V
Frequency (Hz)	50Hz
Max Power (W)	101W
Max Current Rating (A)	0.78A
Energy Consumption (kWh/24h)	0.59kWh/24h
Energy Consumption AE (kWh)	182.5kWh
Heat Output Estimated at 22°C (kJ)	2124kJ
Adjustable Feet (Number & position)	4 All adjustable

Internal Specification	
Number of Shelves	4
Number of Shelf Positions	5
Shelf Size 1 (mm) (WxDxH)	496mm x 405mm x 6mm
Shelf Size 2 (mm) (WxDxH)	493mm x 265mm x 6mm
Shelf Material	White Polythene Coated Wire
Number of Drawers	n/a
Number of Drawer Positions	n/a
Drawer Size 1 (mm) (WxDxH)	n/a
Drawer Size 2 (mm) (WxDxH)	n/a
Drawer Material	n/a
Number of Baskets	1
Basket Material	Polythene Coated Wire
Basket Size (mm) (WxDxH)	485mm x 255mm x 171mm

Weights & Dimensions	
Unit External Width (mm)	595mm
Unit External Depth (mm)	600mm
Unit External Height (mm)	850mm
Unit External Depth Inc Handle (mm)	660mm
Internal Dimensions (mm) (WxDxH)	495mm x 480mm x 720mm
Product Weight (kg)	34.79kg
Packed Width (mm)	640mm
Packed Depth (mm)	660mm
Packed Height (mm)	910mm
Packed Weight (kg)	39.10kg
Country of Origin	China



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.