

Laboratory Plus - ATEX Approved Bluetooth Enabled Refrigerator

394L Laboratory Plus Upright Fridge with Solid Door

Product code LSF400BT-ATEX

SKU code 444411365



ATEX Zone 2, Category 3 Certified (Internal Cavity)

The **internal cavity** of this unit is certified in accordance with ATEX Directive 2014/34/EU for use in **Zone 2, Category 3** environments. This means it is suitable for areas where explosive atmospheres are unlikely during normal operation but may occur briefly and infrequently. All **external components** are classified as being in a non-hazardous (safe) area and therefore do not require ATEX compliance.

Key features


- ✓ Temperature range: of +2°C to +10°C
- ✓ Air temperature probe

- ✓ Audible and visual temperature alarms
- ✓ Factory fitted 16mm porthole for secondary temperature monitoring probes



Bluetooth Connected Controller

Our upgraded controller is Bluetooth enabled and allows for simple, contactless monitoring of temperature and alarm history, via the official Lec Medical app.

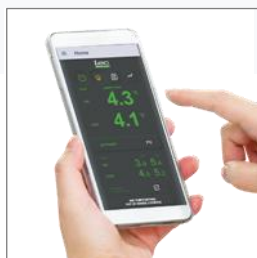

Lec Medical Smartphone App

Contactless local monitoring of air temperature via our bespoke smart phone app. Downloadable temperature reports and alarm history.



Fitted Lock

Integrated lock, comes with two keys for added security and access control, to keep contents safe.



Spark-free

We have removed items such as fans and lights so nothing in the interior of the appliance has a risk of igniting flammable liquids.



Anti-Microbial Handle

Our handles have an anti-microbial powder coating which inhibits the growth of microbes to help prevent cross contamination.

Specification: 394L Laboratory Plus Upright Fridge with Solid Door

Product Information	
Brand	Lec Medical
Category	Laboratory
Sub-Brand	Laboratory Plus
SKU Code	444411365
EAN Code	5052263113657
GDMA Code	LSFSR400BT-ATEX
Description	Free Standing Upright

Door Specification	
Door Finish	Painted Steel / HIPS*
Door Type	Solid
Door Style	Hinged, Reversible
Door Liner	Balcony
Number of Door Trays	7
Number of Door Tray Positions	13
Door Open Depth (mm)	1150mm
Replacement Door Seal	Yes

Product Features	
Gross Capacity (litres)	394L
Net Capacity (litres)	394L
Temperature Range	+2°C to +10°C
Open Door Sensor	Yes
Key Lock	Yes, 2 Keys
Digital Lock	No
Porthole	Yes
Porthole Size (Ø Diameter)	16mm
Porthole Location	Rear Bottom (Right)
Controller Type	Digital BT, External
Temperature Display on Panel	Yes
Temperature Alarm	Yes
Remote Alarm Contacts	Yes
Audible Alarms	Yes
Visual Alarms	Yes
Interior Light	No
Light Type	n/a
Light Activation	n/a
Fan - Type	No
Spark Free	Yes
Defrost setting	Auto (Off Cycle)
Battery Back Up/Hours	No
Bluetooth	Yes
NFC	Yes
WIFI	No
GSM-4G	No
Number of Probes	1
Probe Type	Air
Anti Microbial Handle	Yes

Technical Information	
Noise Level (dB)	36.5dB
Colour	White/green
External Finish	Painted Steel
Internal Finish	HIPS*
Type of plastic	SRL800 (HIPS/ABS mix)
Ambient Temperature Range	16°C to 32°C
Refrigerant Type / Weight (g)	R600a/38g
CFC Free	Yes
HFC Free	Yes
HCFC Free	Yes
Ammonia Free	Yes
Insulation Thickness (mm)	45mm
Insulation Type	High Density Polyurethane Foam
Blowing Agent	Cyclopentane
Certification	CE/UKCA
ATEX Zone Classification	Internal Cavity - Zone 2, Cat 3

Power/Installation	
UK 3 Pin Plug Power Cord	Yes
EU Schuko 2 pin plug Power Cord	Yes
AU Type I Plug Power Cord	No
Voltage Supply (V)	220-240V
Frequency (Hz)	50Hz
Max Power (W)	66W
Max Current Rating (A)	0.49A
Energy Consumption (kWh/24h)	0.58kWh/24h
Energy Consumption AE (kWh)	211.7kWh
Heat Output Estimated at 22°C (kJ)	2088kJ
Adjustable Feet (Number & position)	2 adjustable (front)

Internal Specification	
Number of Shelves	7
Number of Shelf Positions	12
Shelf Size 1 (mm) (WxDxH)	495mm x 325mm x 6mm
Shelf Size 2 (mm) (WxDxH)	487mm x 185mm x 6mm
Shelf Material	White Polythene Coated Wire
Number of Drawers	n/a
Number of Drawer Positions	n/a
Drawer Size 1 (mm) (WxDxH)	n/a
Drawer Size 2 (mm) (WxDxH)	n/a
Drawer Material	n/a
Number of Baskets	1
Basket Material	White Polythene Coated Wire
Basket Size (mm) (WxDxH)	485mm x 180mm x 171mm

Weights & Dimensions	
Unit External Width (mm)	595mm
Unit External Depth (mm)	600mm
Unit External Height (mm)	1925mm
Unit External Depth Inc Handle (mm)	660mm
Internal Dimensions (mm) (WxDxH)	500mm x 480mm x 1775mm
Product Weight (kg)	63.63kg
Packed Width (mm)	640mm
Packed Depth (mm)	635mm
Packed Height (mm)	1995mm
Packed Weight (kg)	70.4kg
Country of Origin	China



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.